ZIRKONUS

Implantatsysteme

SURGICAL MANUAL

ZIRKONUS

IMPLANTATSYSTEME

ZIRKONUS Implantatsysteme GmbH & Co. KG
Bahnhofstrasse 18 – D - 71034 Böblingen
GERMANY

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

Content

1.	Fundamental Principles of the System	5
1.1.	Geometric structure	5
1.1.1.	System components of the multi-part system	6
1.1.2.	System components of the one-piece system	9
1.1.3.	Replacement of multi-root molars with OVERSIZE implants	11
1.2.	Supraosseous principle – placement of the implants in the bone – comparison of the systems	13
1.3.	Placement of implants in the jaw	16
1.4.	Intended use	17
1.4.1.	Purpose	17
1.4.1.1.	Implantatsystem	17
1.4.1.2.	Abutments	17
1.4.1.3.	Cap	17
1.4.1.4.	Gingival abutments and healing caps	17
1.4.1.5.	Planungsbasiselements	17
1.4.2.	Expected clinical benefit	17
1.4.3.	Indication	17
1.4.4.	Contraindication	19
1.4.4.1.	Absolute contraindication	19
1.4.4.2.	Relative contraindication	19
1.4.4.3.	Local risk factors (= relative local contraindications)	19
1.4.4.4.	Temporary contraindication	19
1.4.4.5.	Implant-dependent contraindications	19
1.4.4.6.	Product-specific contraindication ZIRKONUS abutments	20
1.4.4.7.	Product-specific contraindication caps	20
1.4.4.8.	Product-specific contraindication gingival abutments and healing caps	20
1.4.4.9.	Product-specific contraindication of planning basic elements	20
1.4.5.	Risk Factors / Warnings	20
1.4.5.1.	The risk factor	20
1.4.6.	Possible complications	20
1.4.6.1.	Possible procedural complications	20
1.4.6.2.	Implant-related complications	21
1.4.7.	Implant-dependent warnings	21
1.5.	Application, availability, precautions and documentation	22
1.6.	Quality, Warranty and Liability and Development	22

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

2.	The ZIRKONUS Implant System/System Components/System Description	22
2.1.	Materials used	22
2.2.	Structure of the implant/abutment/dental crown/bridge construction	22
2.3.	System Components.....	23
2.3.1.	Components for sterilization.....	23
2.3.1.1.	Tray with inserts (in planning).....	24
2.3.1.2.	Film shrink-wrapping method	26
2.3.2.	Instruments for dissection.....	28
2.3.3.	Instruments for insertion.....	31
2.3.4.	Instruments / Parts for Prosthetics	32
2.3.4.1.	Caps / Healing caps	32
2.3.4.2.	Screwdriver for ceramic screw	32
2.3.4.3.	Rotary wrenches for steel screws.....	33
2.3.4.4.	Screw	34
2.3.4.5.	Planungsbasiselemente (PBE)	35
2.3.4.6.	Anchor Screw-in and Screw-Out Tool	36
2.3.4.7.	Planning template	37
2.4.	Color coding.....	38
3.	Procedure of the treatment.....	38
3.1.	Bone preparation	38
3.2.	Detachment from the packaging / attachment of the tool in the case of one-piece and multi-part systems.....	40
3.3.	Insertion of the implant into the prepared bone cavity.....	41
3.4.	Exposing.....	41
3.5.	Support / fastening of the temporary supply:	41
3.6.	Model production / procedure of the production of a superstructure.....	42
3.7.	Insertion of the abutments / scaffold construction, replacement of the anchor part.....	44
3.8.	Insertion of the sleeve/abutment/temporary or superstructure	45
4.	Reprocessing of surgical instruments and system components	46
5.	Transport and storage of implants	46
6.	Marking / Symbols	47
7.	Specific error possibilities when dealing with the system	48
8.	Imprint	50
9.	Grounds:	51
9.1.	Overview Implant Variants and Sizes	51
9.2.	Overview of multi-part implants	52

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

9.3.	Overview Implants One-Piece	53
9.4.	Overview of rotating instruments.....	54
9.5.	Overview of depth gauges and thread cutters	55
9.6.	Overview Cap - Gingival Abutment	56
9.7.	Overview of planning basic elements	57
9.8.	Overview of impression posts.....	58
9.9.	Overview of impression posts on implants	59
9.10.	Overview of Manipulating Implants	60
9.11.	Insertion process.....	61
9.12.	Overview Scanbody	62
9.13.	Overview Scanbody on Manipulating Implant	63
9.14.	Overview of System Components.....	64
9.15.	A 200 - Instructions for use and preparation.....	65
9.16.	A 201 - Instructions for use implant system	67

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

1. Fundamental Principles of the System

1.1. Geometric structure

The dental ZIRKONUS implant system is a metal-free, one- and multi-part implant system. The coupling geometry between implant, abutment and crown is patented in the multi-part system. It is a ceramic-friendly design in which the dimensioning of the component geometries is optimally adapted to the material properties of the ceramic.

The one-piece and multi-part system differs in that in the one-piece system, the abutment is fused with the implant body to form a complex.

The surgeon who is dealing with the ZIRKONUS implant system for the first time must therefore be aware that there are fundamentally different principles or procedures due to the completely different geometric structure compared to a metal construction.

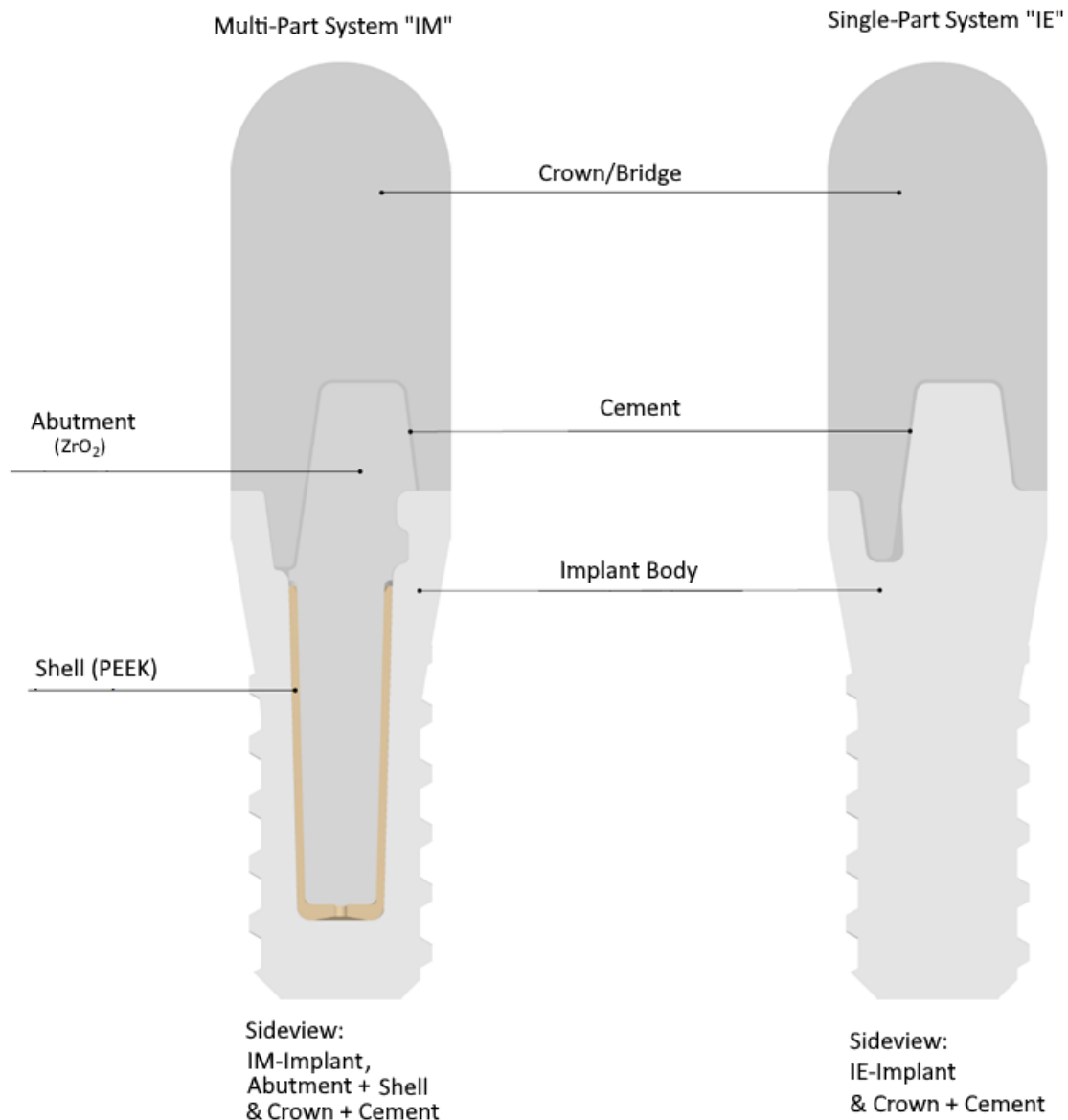
Fig. 1.1-01

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

1.1.1. System components of the multi-part system

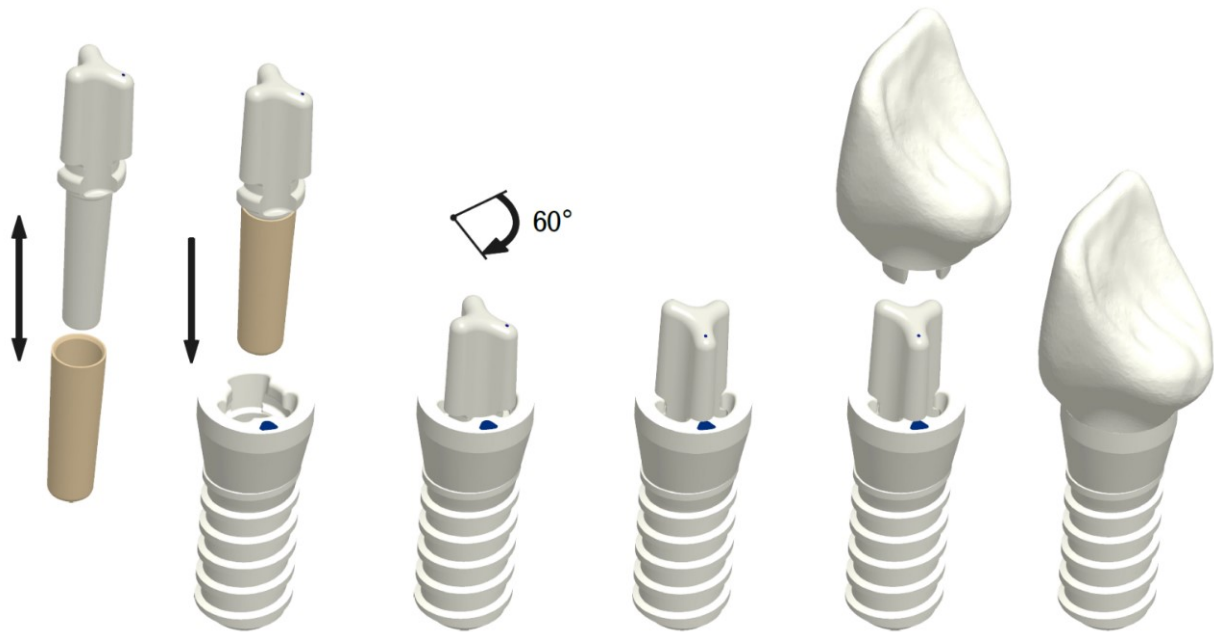
Fig. 1.1-02

According to EU/US patent (EU No. 1992304/US No. 7,726,969 B2), the implant body and the abutment are coupled by a Renk connection (rotation of the abutment by 60 degrees when rotated to the right).

When screwed in, the abutment lowers to the end position at 5°. The rotation is then secured by placing the tooth crown.

The abutment is made of zirconium dioxide, and a PEEK sleeve is required to transmit force from the cone to the inner wall of the implant.

No screw channels are necessary for the superstructure, not even for the abutment. This allows much more delicate abutments and scaffolding constructions to be realized. Cavities in the implant-abutment-crown complex are also minimized (similar to a 3D puzzle), so that there is no exciter reservoir as in metal constructions.

All secondary parts or abutments can be changed as often as desired in the ZIRKONUS implant system. The crowns are cemented onto the implant/abutment complex. Only after the tooth crown has been removed can a reverse rotation with removal of the abutment take place.

There is a color marking on the shoulder of the implant. It enables the precise insertion and replacement of secondary parts. The colour indicates the diameter of the implant (see 2.4.). All instruments for the same implant diameter are equally color-coded.

The instruments are designed in such a way that ceramic surfaces cannot come into contact with metal. This ensures that they are not contaminated by metal abrasion.

The implant diameters \varnothing 4mm (REGULAR) and \varnothing 5mm (WIDE) as well as \varnothing 6mm (LARGE) are pre-cut in the jawbone, while the implant diameters (core cylinders) \varnothing 3.2mm (SMALL), \varnothing 7mm (OVERSIZE), \varnothing 8mm (OVERSIZE X), \varnothing 9mm (OVERSIZE 2X) and \varnothing 10mm (OVERSIZE 3X) are self-tapping implants.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

Ceramic caps and ceramic gingival abutments are fixed exclusively with zirconium oxide screws.

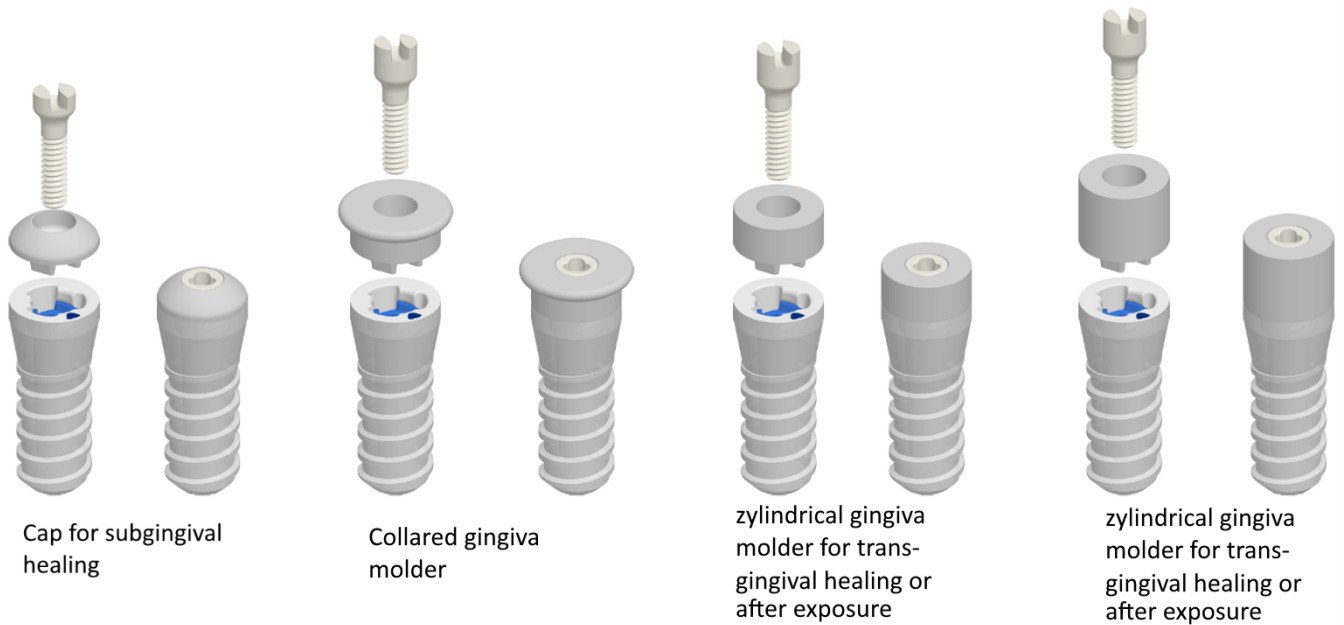
Fig. 1.1-03

The implant is delivered with a pre-assembled anchor part. The anchor part remains in the implant during the healing phase and allows the fixation of secondary parts (closure cap, gingival former, impression posts, planning base elements, etc.).

Fig. 1.1-04

During the healing phase, the anchor part initially remains in the implant body for the fixation of secondary parts. Only at the definitive restoration is the anchor part replaced by the ceramic abutment including PEEK sleeve. (see Fig. 1.1-5)

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

Figure 1.1-05 to Fig. -07 shows the 3D planning in the ZIRKONUS application for a bridge construction in the posterior region. Here, the abutments are designed in height and angle after determining the insertion direction, taking into account the minimum wall thicknesses of the implant crowns.

Fig. 1.1-05

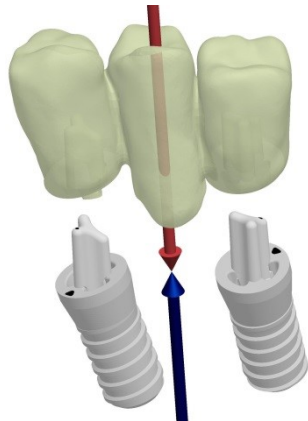
Fig. 1.1-06

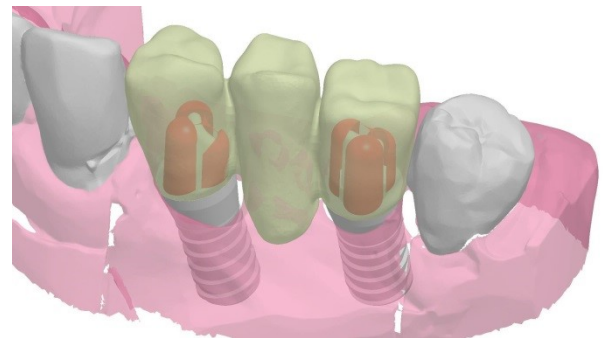
Fig. 1.1-07

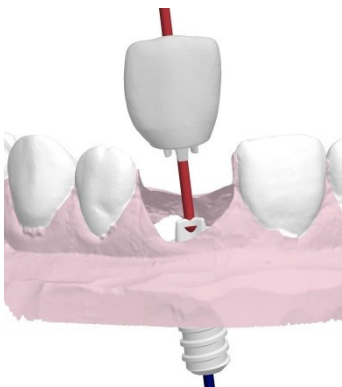
Fig. 1.1-08

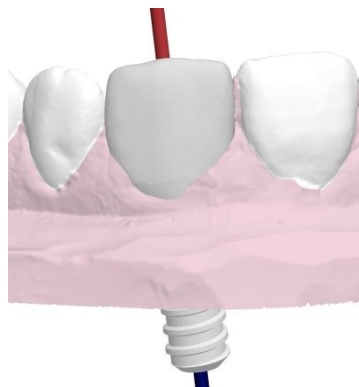
Fig. 1.1-09

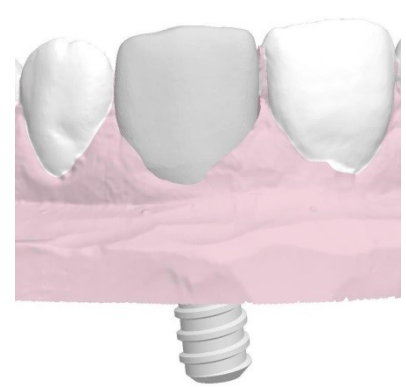
Fig. 1.1-10

Figures Figs. 1.1.08 to Figs. 1.1-10 show the situation of a single crown restoration. Here, too, the insertion direction is important so that the geometry of the adjacent teeth and the antagonists can be taken into account.

The ZIRKONUS implant system is designed in such a way that implant bodies, abutments and scaffold constructions can only be manufactured using the ZIRKONUS CAD/CAM process. Production in any other way is not envisaged and also not possible.

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

1.1.2. System components of the one-piece system

In the case of a single tooth gap, where it is ensured that a secure splint can be attached between the adjacent teeth to prevent unwanted chewing forces during the healing phase, where there has also been no significant vertical bone resorption, so that the intermaxillary distance to the opposing jaw is not too large ($\leq 2\text{mm}$ HDK applied from the upper edge), there are no undercuts due to a problematic insertion direction during the subsequent crown restoration, the practitioner can certainly think about whether to place a one-piece implant instead of a multi-part single implant. In cases that cannot be clearly assessed, however, preference should always be given to the IM system.


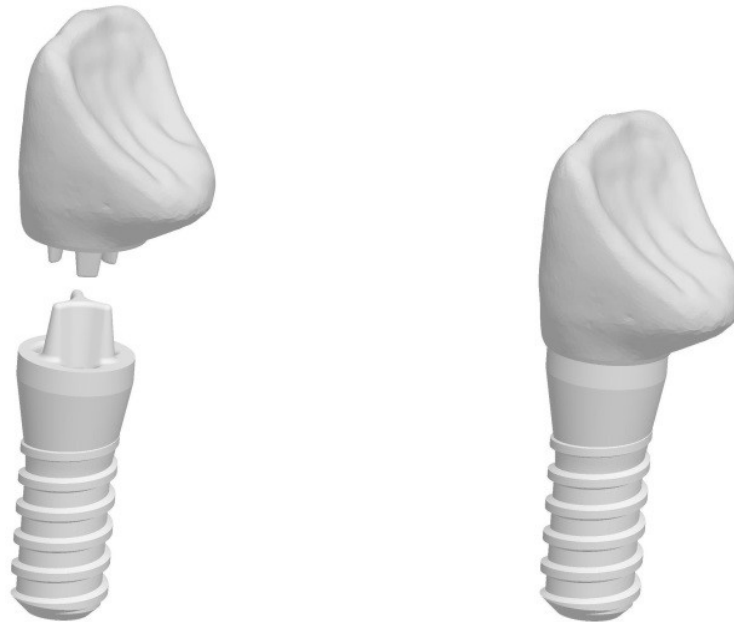
Fig. 1.1-11 Single tooth crown on IE implant

The abutment geometry of the IE implant is fused with the implant body, the crown only needs to be cemented. This is a comparatively inexpensive design for single dentures.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

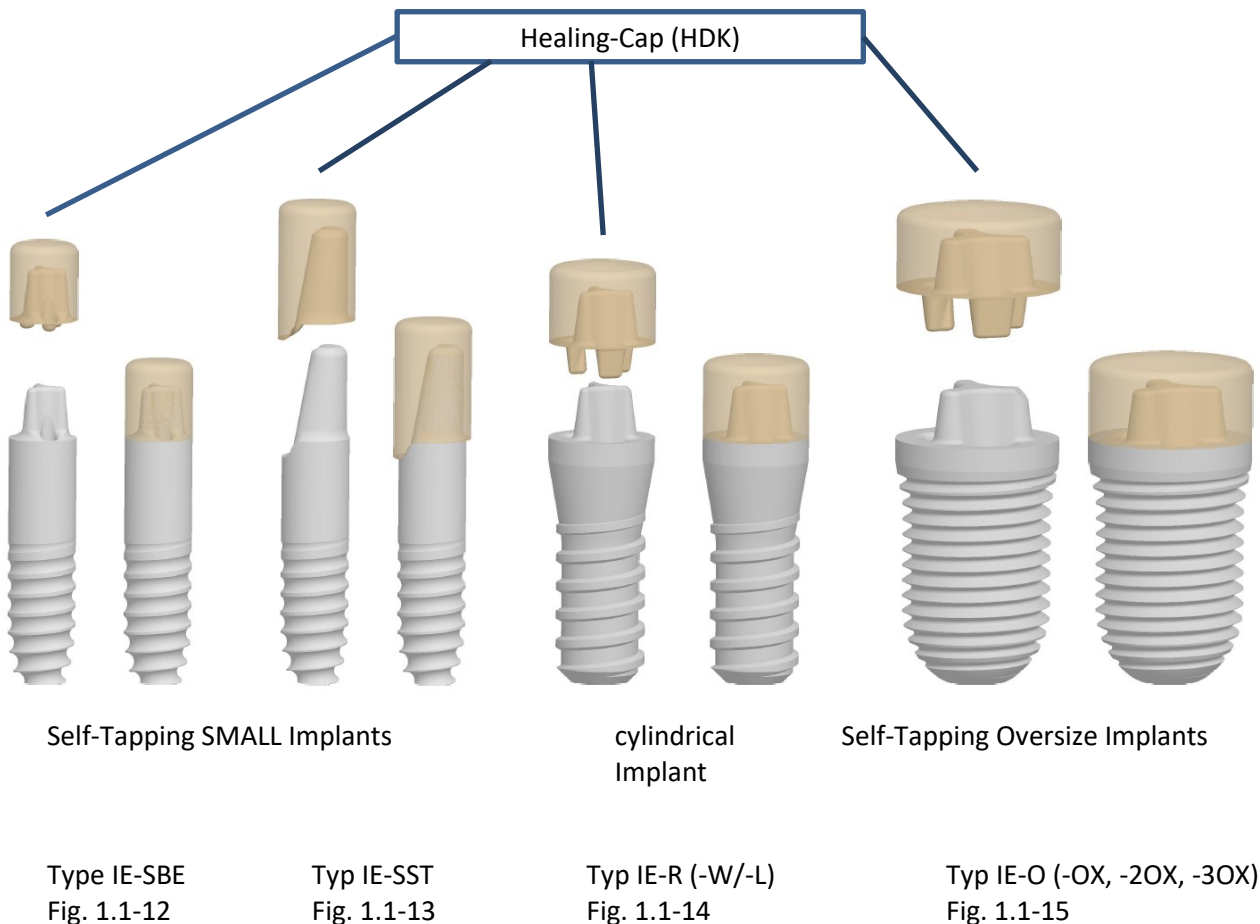The healing cap is placed on the IE implant immediately after insertion into the bone to prevent the gingiva from growing onto the implant shoulder area.

The one-piece system has essentially the same external geometry as the multi-part system.

If stress-free healing is not ensured during the healing phase of the implant due to the situation in the patient's jaw, one-piece implants may not be used.

SMALL implants are available in two different abutment designs (each in two different lengths: L11 (mm) and L14 (mm)).

On the one hand, there is the IE-SBE type (see Fig. 1.1-12), and on the other hand, the IE-SST type (see Figs. 1.1-13).

The IE-SBE type is intended for bridge extension, only in rare individual cases it makes sense as a single tooth prosthesis for a narrow gap in the molar area.

The IE-SST type is used for single tooth restoration in the upper and lower jaw anterior region in narrow gaps, ideally in the case of tooth replacement 12/22, e.g. in the case of aplasia or in the case of replacement of the lower jaw incisors 31/32/41/42.

With this type, it is essential to ensure that the one-sided throat is always directed towards the buccal. It is designed in such a way that the acting chewing force always acts from the back, which has a cylindrical contour parallel to the axis.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

The impression taken after the implants have healed is carried out by the dentist in a conventional way using an individual spoon and plastic impression material using ZIRKONUS impression posts. Unfortunately, enoral scanners are not suitable in terms of accuracy at the moment.

After a plaster model has been produced, it is scanned in the laboratory by a suitable scanner, preferably from the manufacturers: Imetric, 3Shape (900 series), ZfX Evolution. Suitable planning software is Exocad.

After 3D planning, the scan and the desired external geometry are transmitted to ZIRKONUS implant systems via the Internet. There, the abutment and superstructure are produced without intermediate impressions.

1.1.3. Replacement of multi-root molars with OVERSIZE implants

In the shoulder area, the largest ZIRKONUS implant (OVERSIZE 3X) measures 11.5 mm in diameter on the shoulder platform

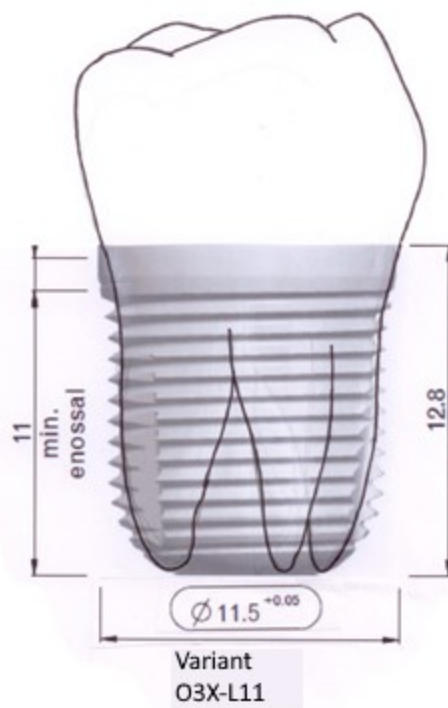
Fig. 1.4-01

In addition to the most commonly available implant diameters of approx. 3.5 to 6.5 mm, ZIRKONUS Implantatsysteme offers ceramic implants with a shoulder diameter of up to 11.5 mm. The corresponding core hole diameter here is a maximum of 9.9mm.

The reason for this is that the anatomical relationships between large and small individuals, men and women, old and young people, and finally also large differences between alveolar processes in the front, premolar and molar regions. In this respect, it is generally surprising that the diameters available on the market since then are usually only between 3.5 and 6 mm.

The so-called oversize implants, available in 4 diameter variants, are specially configured for very cancellous bone, especially for the maxillary molar area. In this case, the bone is very cancellous, i.e. it has large parts of bone marrow and little spongiosa, especially in edentulous areas with reduced height. If implants are placed here, they should have a platform at bone level from the outset that corresponds to the lost tooth with a large diameter, so that a harmonious transition from the implant platform to the attached crown is achieved.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

Much more important, however, is that the threads are supported by the lateral compacta layer and not like a screw with a smaller diameter in the soft, flexible cancellous layer, and thus only an often unsatisfactory primary stability can be achieved. In addition, with the reduced bone height often found there with the much larger bone contact area, enormous primary stability can still be achieved with a larger circumference.

During the definite chewing load in the healed state of the implant, the surface pressure on the bone is reduced and thus reaches values in the order of a multi-root tooth.

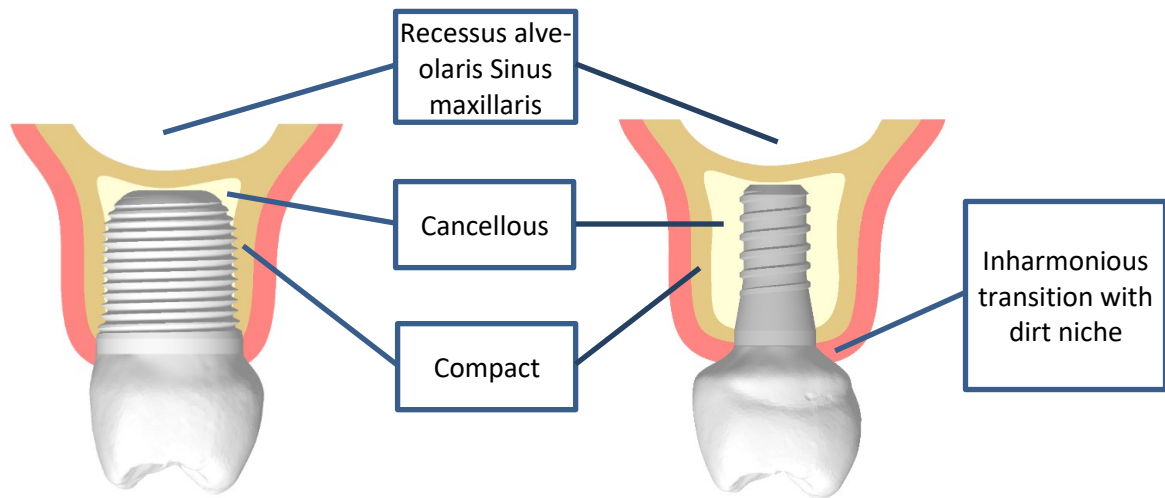
Fig. 1.4-02: Cross-sectional images in the upper chival area with inserted oversize implant and oversize implant conventional screw with D=6mm implant

The threads can be supported from the inside on the compacta layer. The shoulder plateau ensures a harmonious transition to the dental crown.

In addition, it must be noted that the highest chewing force is present in the molar area, and therefore there is the highest risk of breakage for implant bodies and implant abutments. However, the more massive the implants are there, the better the forces are distributed over the larger surfaces (shoulder surface and sheath surface), resulting in a reduction in surface pressure and an overall reduction in compressive stresses, making them enormously resilient. A fracture of the implant bodies can be almost ruled out.

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

1.2. Supraosseous principle – placement of the implants in the bone – comparison of the systems

The white implant material can be located directly under the mucous membrane and you don't have to worry that it will shine through the mucous membrane gray-blue, as is the case with titanium.

Therefore, you don't have to hide the material in the bone as with the bone level principle. As a rule, the alveolar process bone should be fully preserved, regardless of its height, without having to mill it down to the wide plateau required for bone-level implants.

Nevertheless, the surgeon must always pay strict attention to the fact that the bone reaches the level of the implant shoulder everywhere just below it and that it is only exposed in places on the side of the alveolar process (see also Fig. 1.2-03). Either he places the implants deep enough, or he performs augmentative measures on a narrow alveolar process to ensure this. At the highest point of the alveolar insertion, the bone protruding above the implant shoulder level should be left to support the mucosa.

Due to the support of the abutment/crown complex on the implant shoulder, the ZIRKONUS implant system does not cause expansion and splitting between the abutment and the implant body as is the case with thin-walled bone-level implants.

With ZIRKONUS implants, the bone can therefore easily be designed to run pointedly to the implant shoulder; there is no need to fear that the implant wall of ZIRKONUS implants will be stretched and damage the bone.

In reality, it can often be observed that in the first few months after the implants are placed, there is a vertical bone loss of 0.5 – 1 mm. This is irrelevant in the case of the tissue-friendly and white zirconium dioxide (ZrO₂), as it is often observed that the mucous membrane adheres directly to the ceramic and there is no discoloration of the marginal gingiva. Bone resorption of up to 3 mm in relation to the implant shoulder level is also not a problem for the stability of the implant system. This has been proven with dynamic load tests in accordance with DIN EN ISO 14801. There are no undercuts as with bone-level implants with dirt niches with ZIRKONUS implants, but the platform switching principle is still practiced because the manufacturing tolerances of the implant diameters are in the plus and the contact surfaces of the implant crowns are in the negative.

The ZIRKONUS implant system thus takes an extreme opposite position to the bone level principle. Nevertheless, the surgeon must of course keep in mind that the implant is ALWAYS marginally covered by mucous membrane at the lowest point of the bone level and will remain so later! The implant shoulder should be just above bone level if the bone is level. This ensures that there is no separation between the implant body and the implant crown at bone level, and that the cylindrical part of the zirconium implant body should only end above bone level.

In some cases, however, the bone can remain standing above the shoulder level of the implants without any problems. This is advantageous for supporting the inter-implant mucosal papillae.

For orientation during the milling process, the upper edge of the depth marking (black marking, see Fig. 2.3-11) on the cone-shaped moulding cutter shows the practitioner a level above the bone surface of 0.5mm, so that the implant is then placed with the shoulder just above the flat bone level.


In the case of bone-level implants, the author has never made the observation that bone apposition has occurred after the storage of bone substitute material above the gap between the implant shoulder and the abutment.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

The following diagram is intended to show the situation with flat bones as well as with a pointed alveolar process.

Harmonious transition between implant body and dental crown, only one separation, no undercuts.

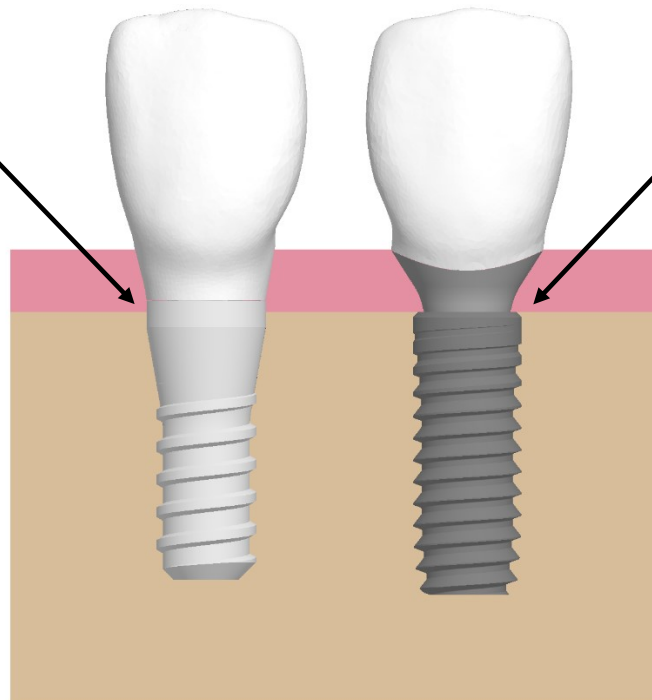Subgingival dirt niche, circularly circumferential volume defect, no mucosal support

ZIRKONUS Implant

Bone Level Implantat

Only one bone level possible, leaving papillary-supporting alveolar process bone standing is not possible. In the case of buccal bone defects, there is a risk of an exposed implant thread.

Fig. 1.2-01 Schematic comparison of mucosal and bone conditions in flat bone

In anatomical situations with flat bone, where total reabsorption of the alveolar proposition has occurred, as is often the case in the mandibular posterior region, there are no serious differences in the placement between ZIRKONUS implants and BoneLevel implants, here the difference is only that the bone is flush with the upper edge of the bone level implant, with ZIRKONUS implants the shoulder surface is just above the bone level. This has the advantage of avoiding irritation of the tissue at the point of passage of the implants on the bone surface, which is present due to material fatigue in the bone level implant principle due to a leaky gap that is often observed.

In the case of narrow alveolar processes, there are fundamental differences between the implant systems when placing the implant bodies. This is shown schematically in the following figure.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

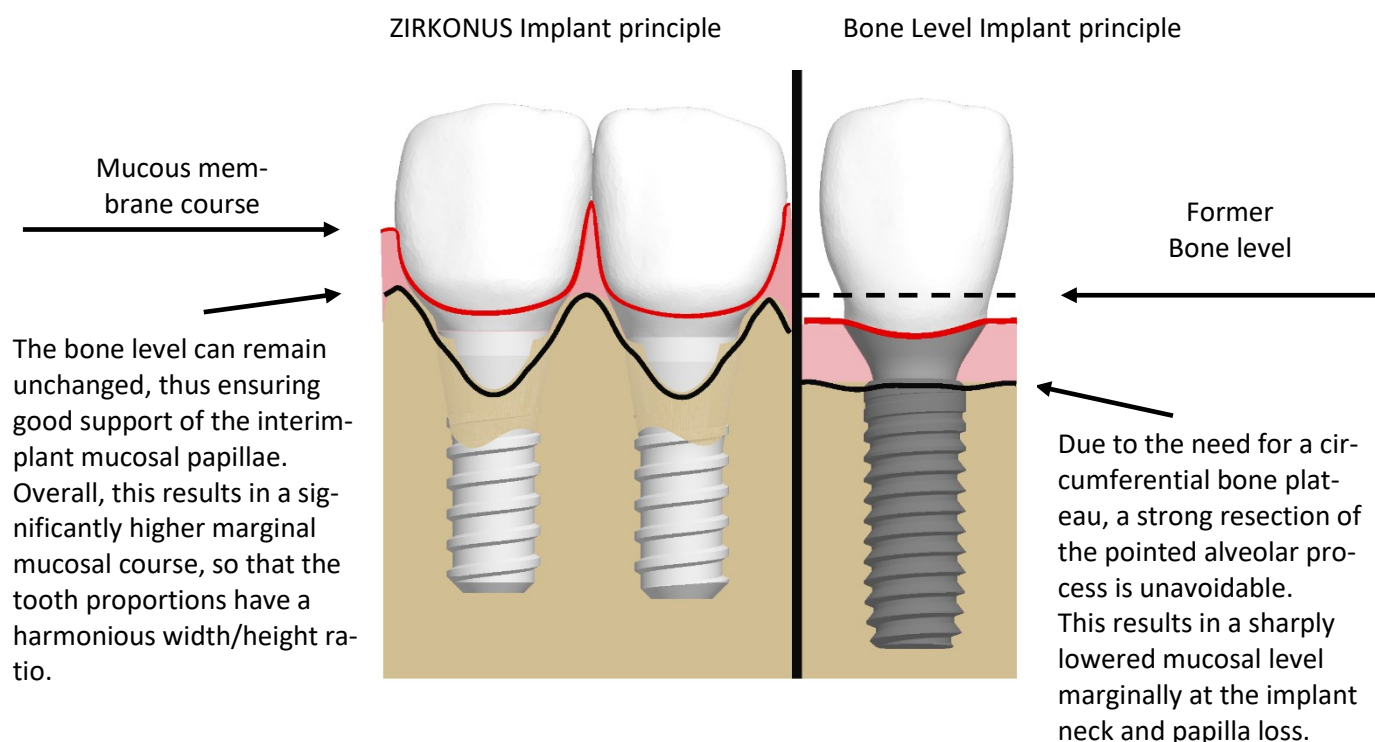
Fig. 1.2-02 Comparison of bone level ZIRKONUS and bone level implant in narrow alveolar process

With the ZIRKONUS implant system, there is no need to mill a plateau on the alveolar process, as is the case with the bone level, in order to ensure that at least 1 mm of bone is still present horizontally around both buccal and orally.

The figure also shows that the bony supporting structures in the bone level implant have been lost as a result, consequently the mucous membrane also sinks, there are no sufficient mucosal papillae, the tooth crown has a significantly longer tooth neck with an unfavorable height-width ratio, which has a cosmetically unfavorable effect. The typical garland-like course of the mucosal papillae can only be realized with the ZIRKONUS implant system. The so-called emergence profile formation in the bone level principle of mucous membrane not supported by the bone can never compensate for this deficiency.

In the case of the bone level implant, two separations are required due to the system: implant / abutment and abutment / crown.

The separation at bone level is, as is well known, extremely problematic because of the gap at bone level, which acts as a bacterial pump. The patient also has no chance of cleaning the circumferential groove at bone level. The second separation in the bone level principle just below the mucosal level inevitably results in an unsatisfactory wide gap, since a cove is not suitable for a ceramic/ceramic connection in terms of manufacturing technology.

Compared to these two weak points of the bone level principle, the ZIRKONUS implant system offers the significant advantage of separation just below gingival level, between implant and crown. Technically, this allows for an almost gap-free transition, because the geometries of the circular radius and shoulder plane are clearly defined, collide at right angles and can therefore be easily implemented in terms of production technology.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

In the ZIRKONUS implant system, the mucous membrane has grown marginally at bone level due to the tissue-friendliness of the material. In the area of the dividing line (implant/crown), which is only just subgingival, it is possible to clean it with dental floss without any problems.

1.3. Placement of implants in the jaw

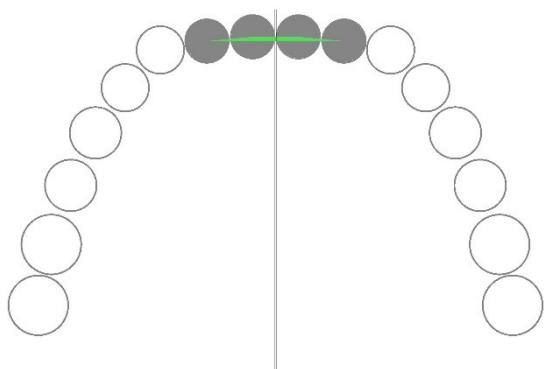
A maximum of one intermediate link should be suspended between 2 implants (distances between the drilled core hole walls may be max. 8mm apart, pendants in bridge constructions may overhang max. 4mm); larger spans are not permitted and are risky due to the risk of fractures or overloading of the implants in the bone. If all teeth are missing in one half of the jaw, it is strongly recommended to place at least 4 implants, preferably 5 or more. In case of reduced bone supply, each tooth should be replaced by an implant, this is especially true in the case of large intermaxillary distance due to the leverage effect. As a general rule, it is important to lock all implant crowns together, this also applies to single crowns of two implants standing directly next to each other.

Blockages with prepared neighboring teeth are possible and useful, especially if bruxism is suspected and unfavorable leverage conditions in the case of long implant crowns.


Only in exceptional situations should larger jaw sections be treated in combination with IE and IM implants, in which case the IM system should always be preferred in cases of doubt. If several IE implants are placed, it is essential to pay attention to parallelism (maximum angular deviation from each other of 15°), especially if a connecting superstructure is provided.

In general, the largest possible implant diameter of the available implants should be chosen for the present alveolar process.

It is important that implants to be treated with blocked crowns are not in a line, but that care is always taken to ensure that a surface is stretched. This is very important in the OC/UK front.

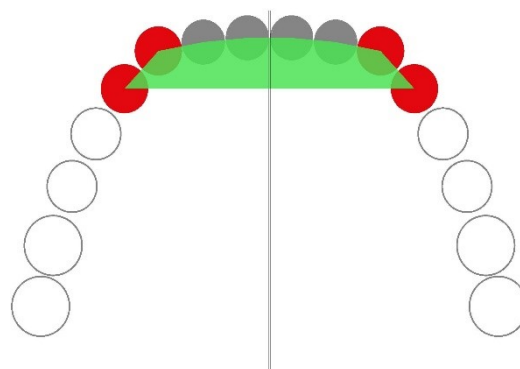If the implants are in line in the front of the upper or lower jaw. In the case of protrusion, there is a lot of stress on the abutments, which poses a risk of fracture

In these cases, it must be urgently advised to include the neighboring pillars in the construction and to block with them.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

1.4. Intended use

1.4.1. Purpose

1.4.1.1. Implantatsystem

The Zirkonus implant system is used for the permanent replacement of the tooth root as well as for the fixation of individual dentures or complete bridge constructions in the event of partial or complete tooth loss.

1.4.1.2. Abutments

ZIRKONUS abutments are components of the IM system (multi-part implants), they serve to ensure that the super-constructions (single crowns, bridges) can be attached to the implants after undisturbed osteointegration. Thus, they are part of the definitive supply.

1.4.1.3. Cap

The closure caps are parts of the IM system, they are suitable for subgingival healing of the implants, if necessary in conjunction with augmentative measures. They are fastened by ceramic screws, which are held in the in situ armature part by the internal thread of the armature.

1.4.1.4. Gingival abutments and healing caps

ZIRKONUS gingival abutments / healing caps are used by temporary fixation on implants surgically inserted into the patient's jawbone to promote a controlled restoration of the soft tissue (gingiva) for later permanent prosthetic restoration by means of a dental crown or bridge.

1.4.1.5. Planungsbasiselements

ZIRKONUS planning base elements are used to support a fixed temporary restoration by means of temporary fixation on implants surgically inserted into the jawbone. The basic prerequisite is resilient and osteointegrated implants of the one-piece IE system and multi-part IM system.

1.4.2. Expected clinical benefit

The expected clinical benefit of the ZIRKONUS Dental Implant System lies in the restoration of chewing function as well as natural oral aesthetics in the case of partial or complete tooth loss.

1.4.3. Indication

The ZIRKONUS Dental Implant System is indicated for fixed dentures (single tooth restoration and bridges of unlimited size), also for immediate implantation.

OVERSIZE implants are primarily suitable for the OK posterior region, in the case of height-reduced and strongly flattened alveolar processes, provided that there is predominantly cancellous bone (type 3 and type 4).

The self-tapping SMALL implants IE-SST and IE-SBE are suitable for narrow tooth gaps with usually narrow alveolar processes.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

The IE-SBE type is intended for bridge extension, only in rare individual cases it makes sense as a single tooth prosthesis for a narrow gap in the molars area.

The IE-SST type is used for single tooth restoration in the upper and lower jaw anterior region in narrow gaps, ideally in the case of tooth replacement 12/22, e.g. in the case of aplasia or in the case of replacement of the lower jaw incisors 31/32/41/42.

Depending on the bone supply (class A and B) and the bone quality (types 1,2,3,4) with different hardness of the bone, the selection of the different types of implants offered must be made by the practitioner. Assistance for the selection of suitable implants for the bone classes and bone quality types can be found in the following table.

<p>Bone quality →</p> <p>Bone quantity ↓</p>	<p>TYPE 1</p> <p>Mainly Compacta, less than 10% cancellous, little elasticity e.g. age-related atrophy or atrophy in the case of non-tooth attachments, problematic when preparing the bone storage due to high hardness</p>	<p>TYPE 2</p> <p>Balanced ratio of Spongiosa to Compacta 30% - maximum 60% Spongiosa good elasticity</p>	<p>TYPE 3</p> <p>Spongiosa predominates 60% - 90%, but still so much compacta that bone can be easily impressed, still elasticity present</p>	<p>TYPE 4</p> <p>Compacta layer is less than 10% in favor of cancellous, compressible, hardly any elasticity no regression after impressing, mostly maxillary posterior region</p>
<p>A</p> <p>Slight to moderate loss of height of the alveolar process, almost complete width of the alveolar process preserved</p>	<p>Preferably Cylindrical Small/ Regular/ Wide/ Large implants possible</p> <p>Rare anatomical conditions</p>	<p>All implant geometries and diameters are possible without restrictions.</p> <p>This type often occurs with recent traumatic tooth loss.</p>	<p>Preferably self-tapping bone-tightening implants of the sizes Oversize, Oversize X, Oversize 2X or Oversize 3X possible, cylindrical implants possible</p> <p>but severely limited indication for small implants</p>	<p>Self-tapping implants of size Oversize, Oversize X, Oversize 2X and Oversize 3X preferred</p> <p>Limited indication for cylindrical implants</p> <p>Contraindication for small implants</p>
<p>B</p> <p>Strong resorption of the alveolar process, at most rudimentary, defects in the jawbone also possible</p>	<p>Prefers cylindrical implants</p> <p>Limited indication for small implants due to risk of fracture during insertion</p> <p>Only for large bone defects Oversize (O, OX, O2X, O3X) implants, otherwise contraindication for OVERSIZE implants</p>	<p>All implant geometries, depending on the desired diameter</p> <p>Good indication for small implants.</p> <p>Only use large bone defects Oversize, Oversize X, Oversize 2X or Oversize 3X implants.</p>	<p>Main indication for self-tapping oversize implants</p> <p>Severely limited indication for cylindrical implants, severely limited indication for small implants due to insufficient primary stability.</p>	<p>Main indication for self-tapping oversize implants</p> <p>Severely limited indication for cylindrical implants,</p> <p>Contraindication for small implants</p>

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

1.4.4. Contraindication

1.4.4.1. Absolute contraindication

- Infected extraction sockets, larger apical osteitis (bone inflammation) and bone inflammation in the jaw area.
- Serious diseases of the bone, metabolism, blood clotting, circulation, heart and immune system, etc.
- Immunosuppression
- Alcohol or drug abuse

1.4.4.2. Relative contraindication

- Radiotherapie
- Schwere Diabetes, speziell juvenile Diabetes
- Hemorrhagic diathesis or anticoagulation
- Nicotinabusus
- OVERSIZE implants: Average compacts in the mandibular posterior region

1.4.4.3. Local risk factors (= relative local contraindications)

- Erosive or bullous diseases of the mucous membrane of the alveolar process
- Bruxism, Parafunctional Habits
- Severe bone resorption

1.4.4.4. Temporary contraindication

- Untreated periodontitis
- Root residues at the implantation site (possibly immediate implantation after extraction)
- Premises Infekt

* The use of the products must be weighed against the risks on a case-by-case basis and carefully considering the expected clinical benefits

1.4.4.5. Implant-dependent contraindications

Due to the usually narrow alveolar process, OVERSIZE implants are not suitable in the anterior and premolar area.

SMALL implants IE-SST and IE-SBE are contraindicated if there is a large tooth gap with a wide alveolar process. The general principle here is that the entire width of the alveolar process should always be used, i.e. that the largest possible implant diameter should always be chosen.

One-piece implants are not suitable for use in conjunction with bone-building measures and are contraindicated for such procedures.

If load-free healing is not ensured during the healing phase of the implant due to the situation in the patient's jaw, one-piece implants may not be used.

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

1.4.4.6. Product-specific contraindication ZIRKONUS abutments

ZIRKONUS abutments may only be used for permanent restoration in conjunction with multi-part implants provided for this purpose by ZIRKONUS Implant Systems and are contraindicated for any use with third-party systems.

1.4.4.7. Product-specific contraindication caps

ZIRKONUS caps are not intended for transgingival healing, gingival abutments and healing caps are available for this purpose, see the following chapter.

1.4.4.8. Product-specific contraindication gingival abutments and healing caps

Gingival abutments and healing caps may only be used in preparation for permanent restoration in conjunction with implants provided for this purpose by ZIRKONUS Implant Systems and are contraindicated for any application with third-party systems.

1.4.4.9. Product-specific contraindication of planning basic elements

ZIRKONUS planning base elements may only be used as a support for a temporary restoration on implants from ZIRKONUS Implant Systems and are contraindicated for any application with third-party systems.

1.4.5. Risk Factors / Warnings

1.4.5.1. The risk factor

- Previous or existing therapy with bisphosphonates (intravenous or oral) for the treatment of metastatic disease, osteoporosis, or osteopenia

1.4.6. Possible complications

1.4.6.1. Possible procedural complications

- Post-operative bleeding
- Trauma to the surrounding anatomical structures
- (Persistent) Ache
- Neuropathies or paresthesia
- Inflammation and/or infection (acute or chronic)
- Fistulas (nasal or maxillary sinus)
- Fracture of the jaw
- Reactive hyperplasia of the gingiva
- Keloid
- Untreated implant, unstable implant, implant loss
- Loss of function
- Progressive bone loss
- Increased exploratory depth
- Persistent peri-implant radiolucency

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

1.4.6.2. Implant-related complications

- Implantatfracture
- Fracture of occlusal materials
- Broken or loosened prosthetic components
- Loss of abutments or screws

1.4.7. Implant-dependent warnings

Classes A3, A4, B1, B3, B4 are critical for self-tapping small implants, and caution is also advised with type 1 bone when inserting bone-densifying conical oversize implants, because a situation could arise whereby the implant cannot be placed in the final position, on the other hand, there are also great difficulties in removing it from the bone because it has seized.

With the implant type IE-SST, it is essential to ensure that the one-sided groove is always directed towards the buccal. It is designed in such a way that the acting chewing force always acts from the back, which has a cylindrical contour parallel to the axis.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

1.5. Application, availability, precautions and documentation

The dental ZIRKONUS implant system is intended exclusively for dentists who are familiar with dental implantology, including diagnosis, preoperative planning, surgical procedure and prosthetic care. Before use, the user must ensure that he has developed and understood the instructions for use, manuals and other product-related information provided by ZIRKONUS. In addition, it is mandatory to attend a training course for the application of the ZIRKONUS dental implant system offered by ZIRKONUS in order to learn the safest technique for the system, since the instructions for use and manuals cannot possibly cover all eventualities in the safe application and replace the personal experience of the tutors.

A product catalogue (attached) is available for the selection of components.

It is the sole responsibility of the user to examine and inform his patients in detail before using ZIRKONUS implants – as usual. In order to reduce liability and ensure optimal implantation success, ZIRKONUS recommends complete clinical, radiological, photographic and statistical documentation of the anamnesis.

1.6. Quality, Warranty and Liability and Development

Development, clinical trials, manufacturing, quality assurance, risk management and market observation are carried out in accordance with Directive 93/42/EEC - Medical Devices Directive and national medical device law.

In principle, there is an exclusion of warranty and liability in the event of improper use of the system components by the user or third parties; this applies in particular to a combination of system components of the ZIRKONUS implant system with third-party products, unless ZIRKONUS has expressly recommended their compatibility. The use of implants and system components of the ZIRKONUS implant system is outside the control of the manufacturer and is the sole responsibility of the user.

The technical advice (oral and written) is carried out in accordance with the state of the art in science and technology when placing the dental ZIRKONUS implant system and its components on the market. It does not relieve the user of the obligation to carefully check their suitability for the intended indications and applications. The procedures for implantation described here are only recommendations, but they are based on many years of experience of the experts of the company ZIRKONUS, from which, however, no assurances or warranty commitments can be derived.

All system components are subject to continuous further development, taking into account the latest progress in the state of the art and new scientific findings. **Changes in construction, design and material are therefore subject to change.**

2. The ZIRKONUS Implant System/System Components/System Description

2.1. Materials used

Only biocompatible materials are used for the system components of the implant system from ZIRKONUS.

For the definitive restoration, zirconium oxide ceramic (ZrO₂), which is approved for surgical implants, and the high-performance plastic PEEK (polyetheretherketone) are used.

In addition to zirconium oxide ceramics, ZIRKONUS instruments are also made of the high-performance plastic PPSU (polyphenylsulfone), surgical, stainless steels and titanium.

(e.g. ratchet / ratchet adapter / elbow shanks / coupling bolts / RA shanks)

2.2. Structure of the implant/abutment/dental crown/bridge construction

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

In the case of implants placed in such a divergent manner, a blocked superstructure can only be realized with the IM system, where the abutments are parallelized for a common insertion direction.

With the IE system, this situation would have to be solved with two single crowns that are not blocked with each other. Blocking of the IE system is only possible up to a total of 15° deviation of the implant axes from each other.

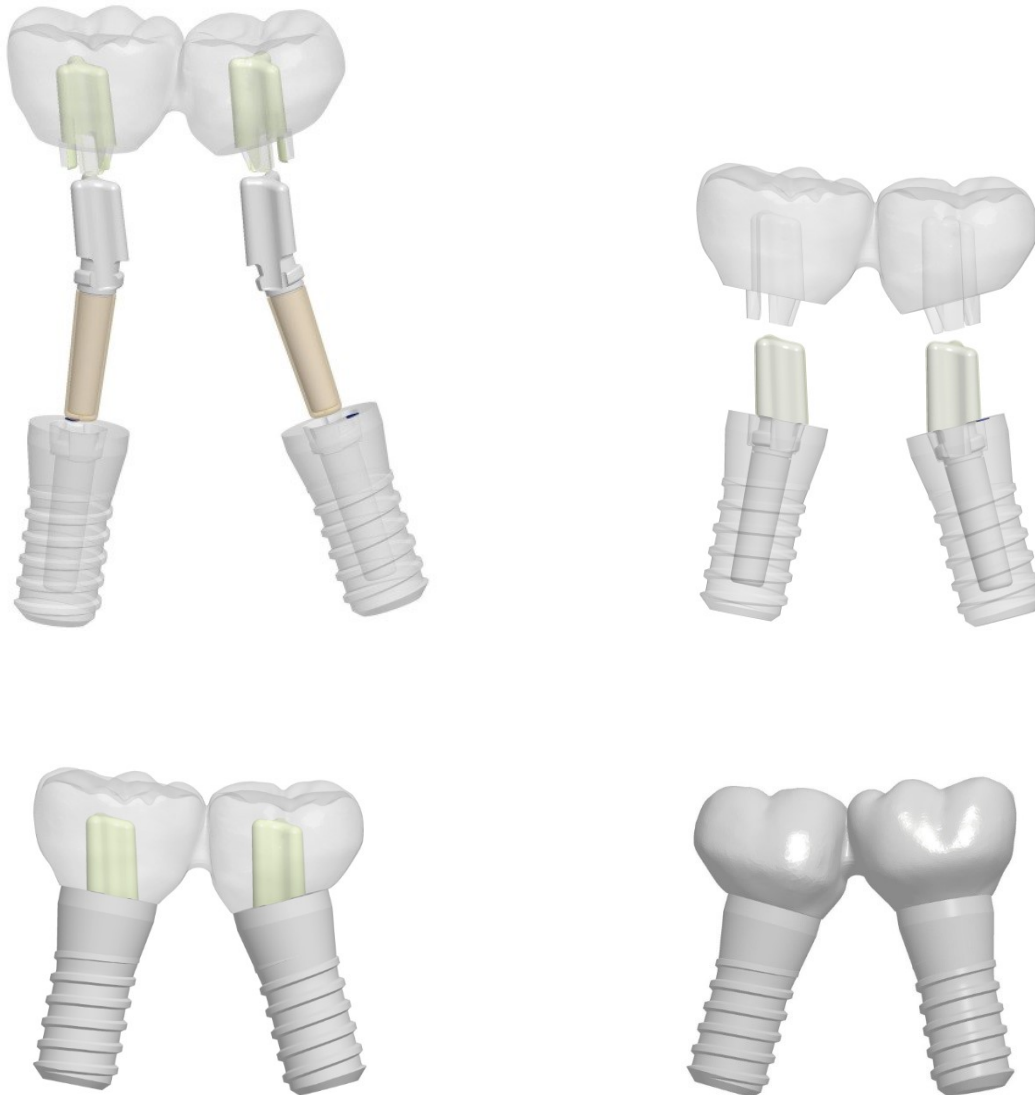
Fig. 2.2-01 Implant, individual abutments, bridge on non-axially parallel implants

A prosthetic restoration with a blocked superstructure is still possible with the IM system up to a divergence of the implant axes to each other of 60°, because the individual abutments can be angled up to 30° in extreme cases, normally axial divergences of 15° per implant should not be exceeded. The height of the abutments can also be individually adjusted according to the situation.

2.3. System Components

2.3.1. Components for sterilization

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

2.3.1.1. Tray with inserts (in planning)

A flat base cassette made of zirconia is provided for steam sterilization of small parts with various removable dividers.

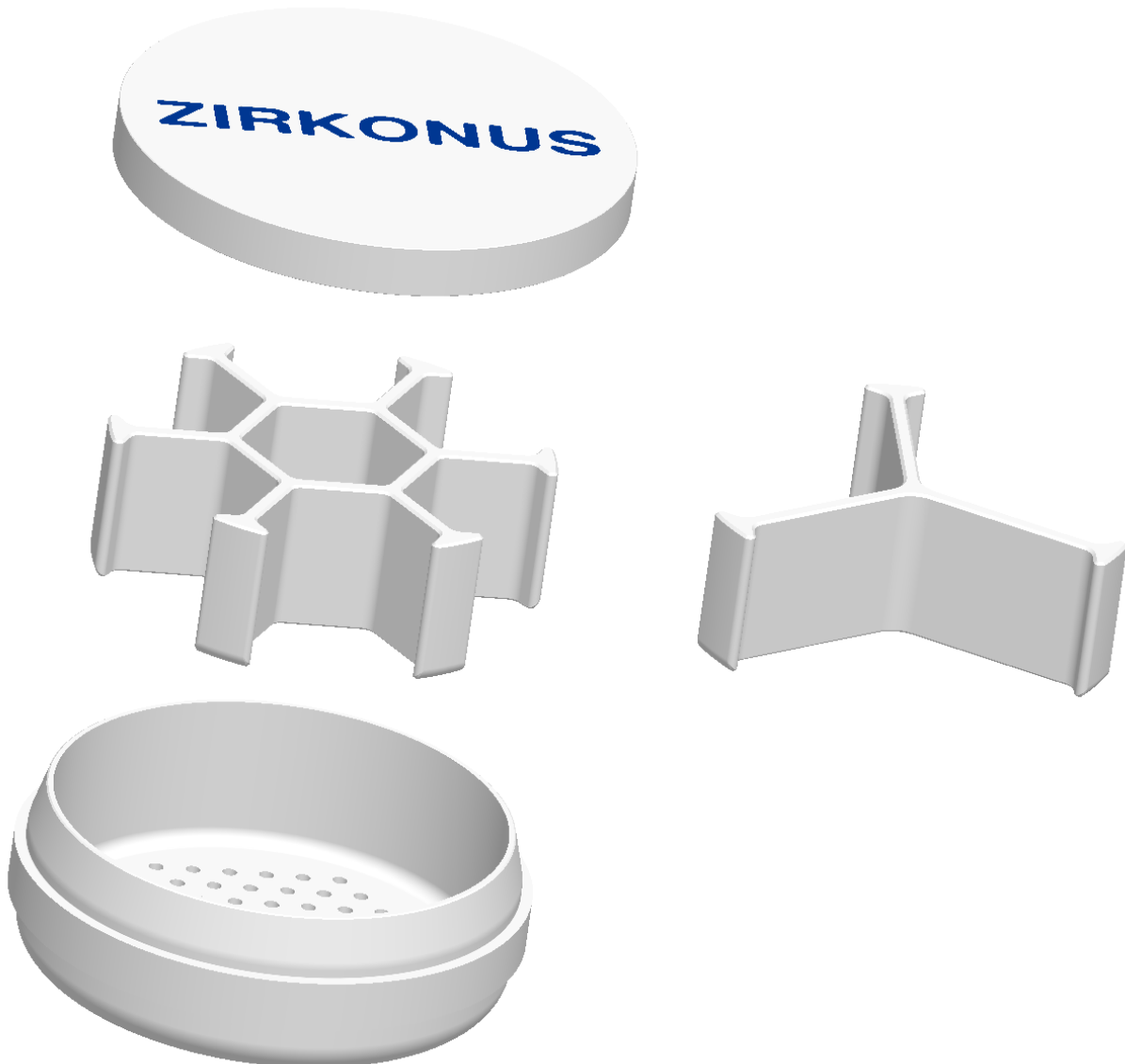
Fig. 2.3-01 Basic cassette made of zirconium oxide with an outer diameter of 80mm, Height 39 mm and inserts

We also provide a cassette with a high lid for instrument carriers. The instrument trays are divided for instruments for small implant diameters: small, regular, wide, large and for large oversize implant diameters.

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

Fig. 2.3-02 Basic cassette made of zirconium oxide with an outer diameter of 80mm and a high lid (total height 48mm) and for instrument carriers:
Left Fig.: Small Implant Diameter Right Fig.: Large Implant Diameter

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

Fig. 2.3-05 2 round ceramic containers (basic cassettes) can be placed next to each other in a sterilization box

The sterilization box should have a minimum interior space of length $L = 165$ mm, width $W = 85$ mm and height $H = 40$ mm. Alternatively, the ceramic containers with their contents can also be shrink-wrapped in sterilization foil.

The inserts in the ceramic basic cassettes are used for the orderly storage, cleaning and sterilization of the ZIRKONUS system components. The colour coding on the insert makes it easier to assign them to the implant variants.

2.3.1.2. Film shrink-wrapping method

We generally recommend shrink-wrapping and sterilizing instruments separately according to materials for the respective operation in sensible combinations.

In the case of drills, care must be taken to ensure that the titanium shanks do not come into contact with the ceramic drill bodies in order to avoid metal abrasion on the ceramic (see Fig. 2.3-06 and Fig. 2.3-07)

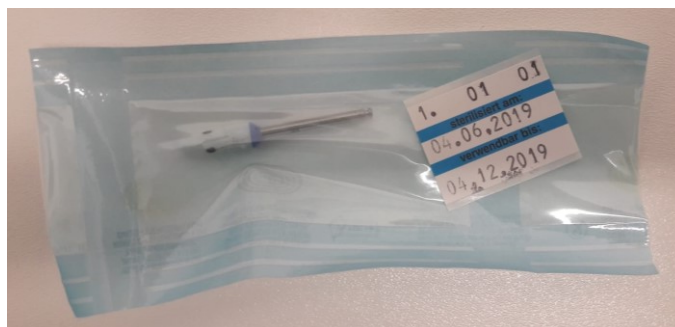
Fig. 2.3-06

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

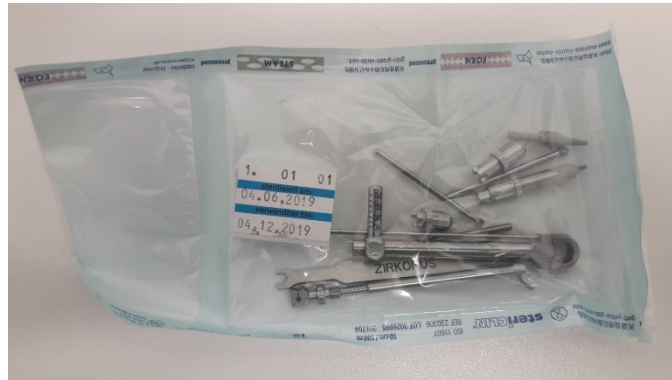
Fig. 2.3-07

When cleaning in the thermal disinfectant, it is also important to ensure that the instruments do not have any contact between the ceramic and the metal. The drills are clamped into the grid of the bracket on the titanium shaft.

Assembled instruments such as thread cutters are dismantled for cleaning and sterilization and only put together before use.

Fig. 2.3-08 Instrument holder for thermal disinfectant

This is ideal for clamping instruments with titanium shafts. This ensures that the ceramic drill bodies have no contact with the metal during processing.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.2. Instruments for dissection

The system components are supplied non-sterile and can be steam sterilized.

The bone drills for cavity preparation are available in the diameter range from 2.5 to 10 mm.

The depth coding allows an orientation with regard to the depth of the bore and an estimation of the bone level in relation to the later location of the implant.

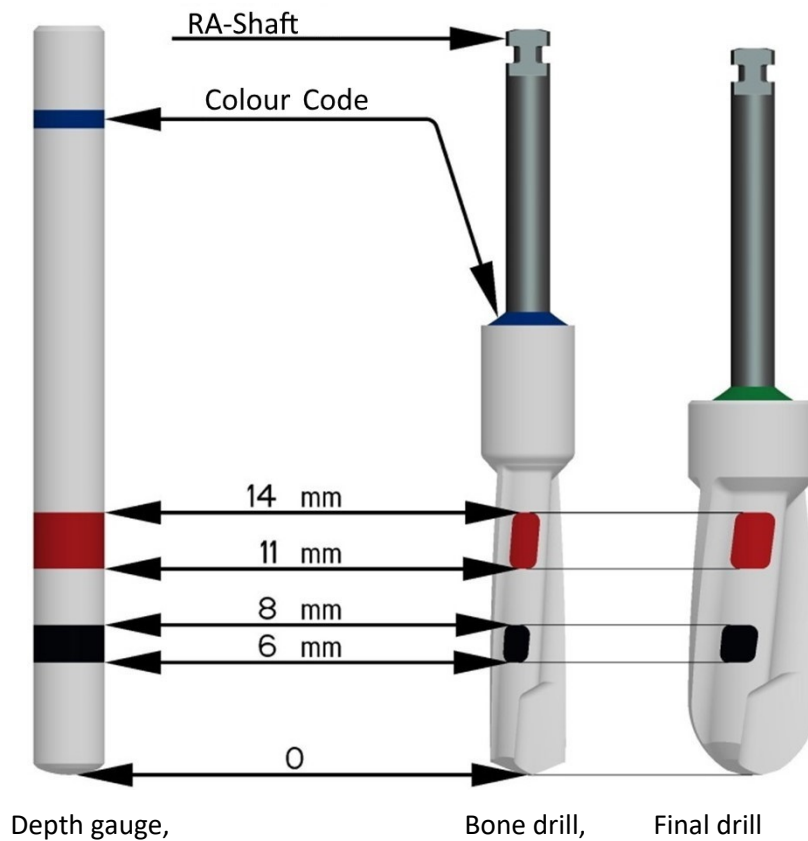
Fig. 2.3-10 Depth gauge,

Bone drill,

Final drill

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

For each desired implant diameter, the pre-drilled hole is then shaped with the moulding cutter. The black marking allows an estimate of the later fit of the implant when milling out with the shape cutter.

Sink the cutter until the bone reaches the lower edge of the mark at its lowest point. The general bone level should reach at least to the top of the mark.

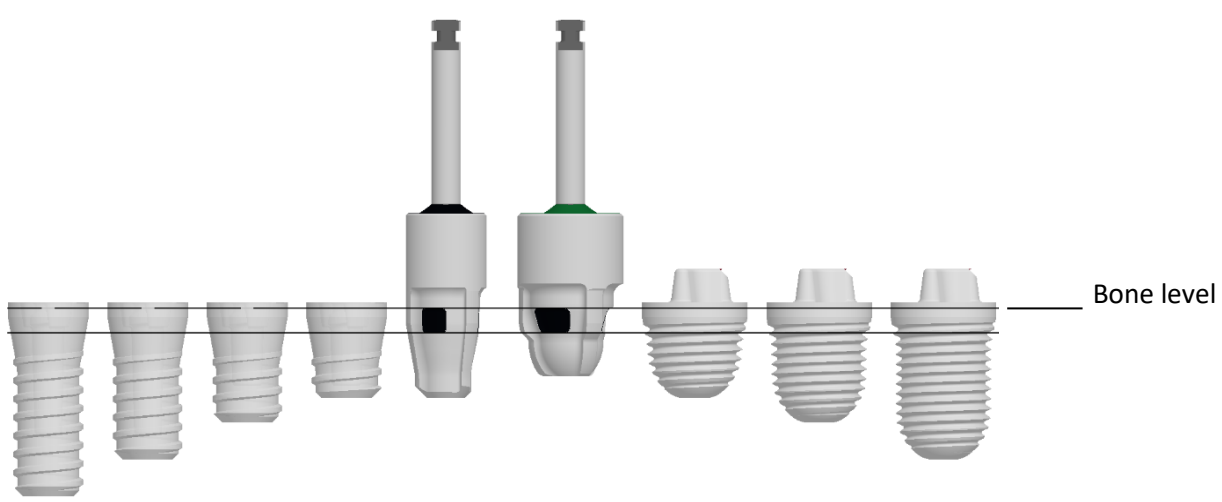
Fig. 2.3.-11: IM-System Moulding cutters IE-System

The instruments for bone preparation are identical for the implant systems.

There are no differences in the external geometry of the implants between the IE and the IM system. In the figure, the IM system is shown on the left as an example for cylindrical implants, on the right the IE system for over-size implants as an example.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

Thread Cutter Components

For the implant variants "Regular, Wide, Large", a thread must be pre-cut in the bone.

The tapping instrument consists of three components.

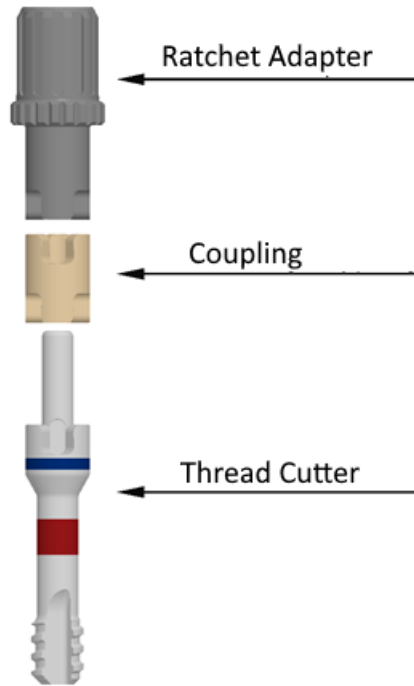
Fig. 2.3-12

The thread cutter (ZrO2) is connected to the ratchet adapter by means of PEEK connection coupling. It can be steam sterilized, and components can be disassembled!

Gewindeschneider						
GS-R	GS-W	GS-L	GS-R with coupling	GS-R with turning tool 03 for ratchet	with turning tool 08	with turning tool 11
						

Fig.2.3-12a

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

2.3.3. Instruments for insertion

After removing the implant holder from the blister primary packaging, the implant is coupled into the insertion tool.

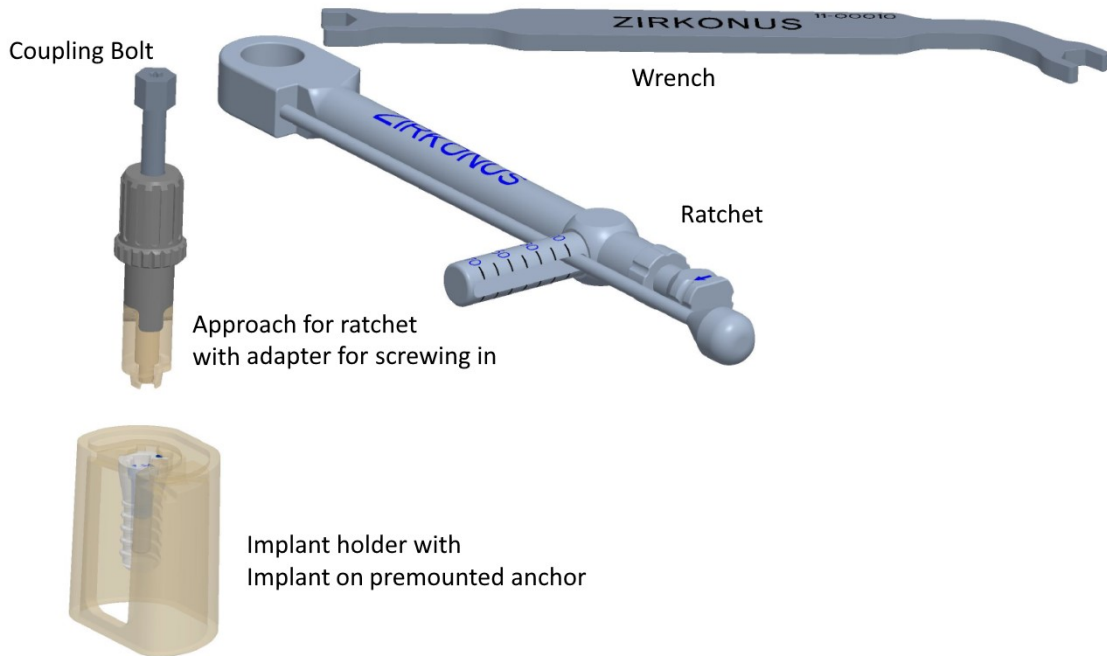
Fig. 2.3-13 multi-part system

Representation of the system components: Screw-in tool consisting of ratchet attachment/screw-in adapter/ratchet/coupling bolt/open-end wrench Implant in implant holder with pre-assembled anchor part

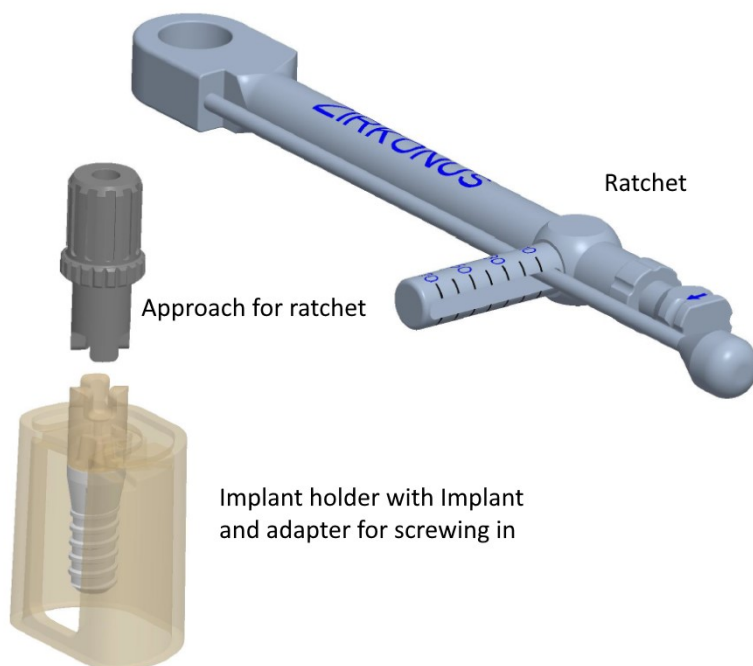
Fig. 2.3-14 one-piece system

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.4. Instruments / Parts for Prosthetics
2.3.4.1. Caps / Healing caps

After the implant has been placed or has been exposed, the following geometries are available for each diameter to close the implant and/or shape the gingiva.

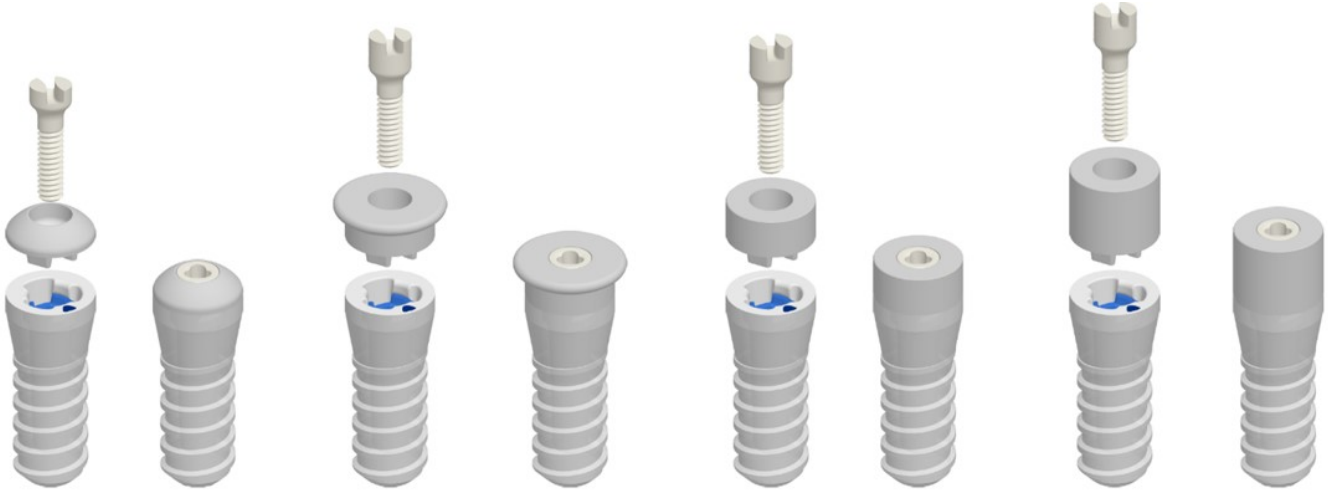
Fig. 2.3-15

Cap for subgingival healing

Fig. 2.3-16

Extended Gingiva Shapes

Fig. 2.3-17

Cylindrical gingival former for transgingival healing or after exposure

Fig. 2.3-18

cylindrical gingival former for transgingival healing or after exposure

2.3.4.2. Screwdriver for ceramic screw

Interchangeable bits made of PEEK are available for the ratchet adapter or the handpiece drive with shank adapter typ 1 (DIN EN ISO 1797-1), screwdriver Bits S06 in different lengths, suitable only for our ceramic screws.

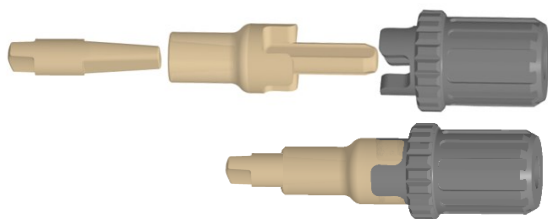
Fig. 2.3-19 a

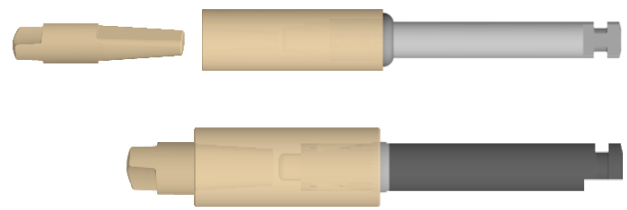
Fig. 2.3-19 b

Screwdriver Bit S06 in different lengths

2,35 x 11mm

2,35 x 15mm

2,35 x 20mm

2,35 x 25mm

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

Fig. 2.3-19 a

Ratchet attachment with swivel wrench adapter for PEEK bits for ceramic screws, Standard screwdriver for ceramic screws.

Fig. 2.3-19 b

Adapter Twist Wrench for PEEK screwdriver Bits with shank adapter typ 1 (DIN EN ISO 1797-1) for coupling to the angle turner tool

Only use PEEK screwdriver for ceramic screw

Fig. 2.3-20 Flat head

Fig. 2.3-21 High head

2.3.4.3. Rotary wrenches for steel screws

For steel screws, the rotary wrench is provided with a hexagonal and connection for the torque ratchet. It is also used for screwing in SMALL implants, see also Fig. 2.2-23.

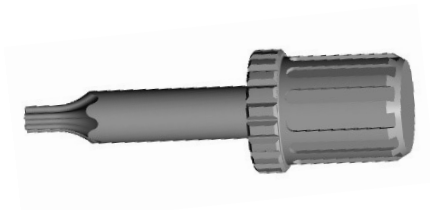
Fig. 2.3-22 Rotary wrench T6

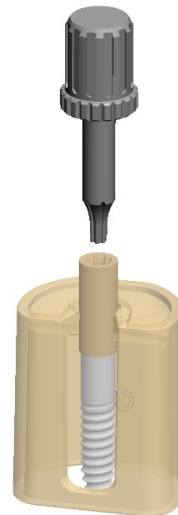
Fig. 2.3-23 Rotary wrench T6 with screw-in coupling part for SMALL implants

Fig. 2.3-24 Contra-angle wrench

The steel wrenches must never be used for ceramic screws.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.4.4. Screw

Fig. 2.3-25 Screw impression post

Fig. 2.3-26 Screw Scanbody

The screw impression post is described in Figure 3.6-01, screw scanbody see Fig. 3.6-02.

Fig. 2.3-27 Screw anchor

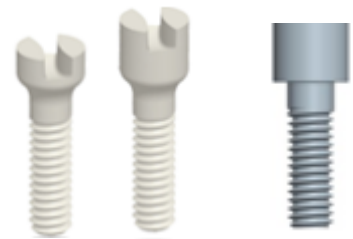
Fig. 2.3-28 Screw PBE (ceramic and steel)

The screw anchor is used in conjunction with the anchor screw-out and screw-in pull-in, see Fig. 3.7-01. It is also used for the screw pull for the IM system, see Figs. 3.2-02 and Figs. 3.3-02.

The Planning Base Element (PBE) screw is used to fix the planning base elements, which are available in different diameters, see Fig. 2.3-30.

Fig. 2.3-29 sleeve remover tool with sleeve short and long

In the event that a definitive abutment inserted with a sleeve has to be removed again, for example in the event of a fracture or during the try-in, and the sleeve does not stick to the cone of the abutment, but is still stuck to the inner wall of the implant, this ceramic instrument (sleeve remover tool) is very helpful to remove the sleeve from the implant, a few turns of the instrument are sufficient, so that the fine thread builds up adhesion in the malleable PEEK sleeve and so that the sleeve can usually be removed easily and quickly.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.4.5. Planungsbasiselemente (PBE)

After the implants have healed, it is possible to support a fixed, temporary restoration on the implants. For this purpose, planning base elements made of PEEK are available, which can be fixed with steel or ceramic screws. It is advisable to place the flattened phase on the side of the strongest tilt in order to achieve an insertion direction without undercuts.

Ceramic screws are available in case there are concerns about metal in the patient concerned. If necessary, metal can also be completely dispensed with during the provisional supply.

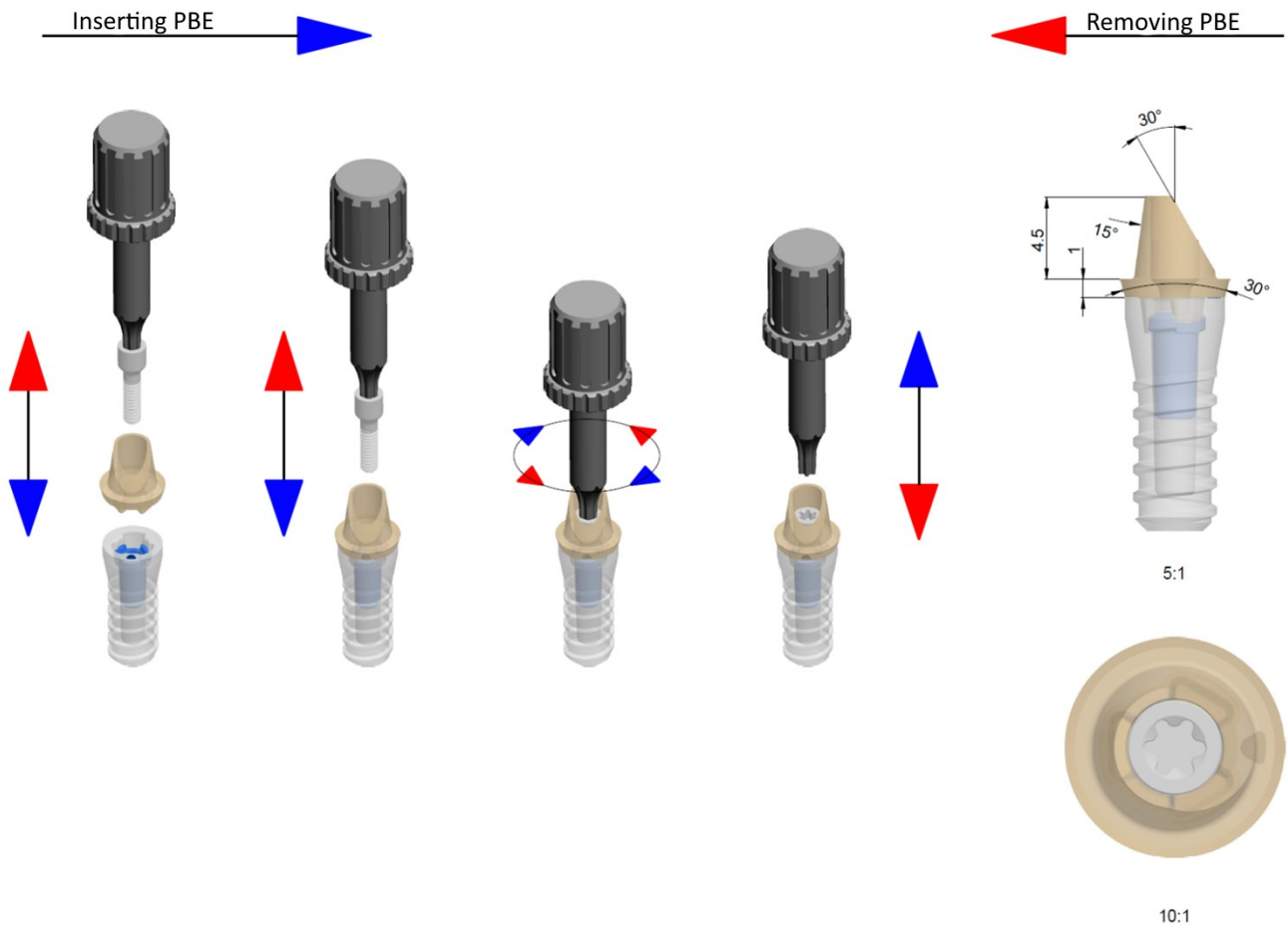
Fig. 2.3-30

The plastic anchor ensures that the steel screw cannot come into contact with the ceramic implant.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.4.6. Anchor Screw-in and Screw-Out Tool

Once all phases of the treatment have been completed, including the impression in the jaw, and the final work on the IM system has been completed with the delivery of the individual abutments and the superstructure but the dental laboratory, the anchor parts must be removed from the implants.

A special instrument, the anchor insertion and removal tool, is available for this purpose.

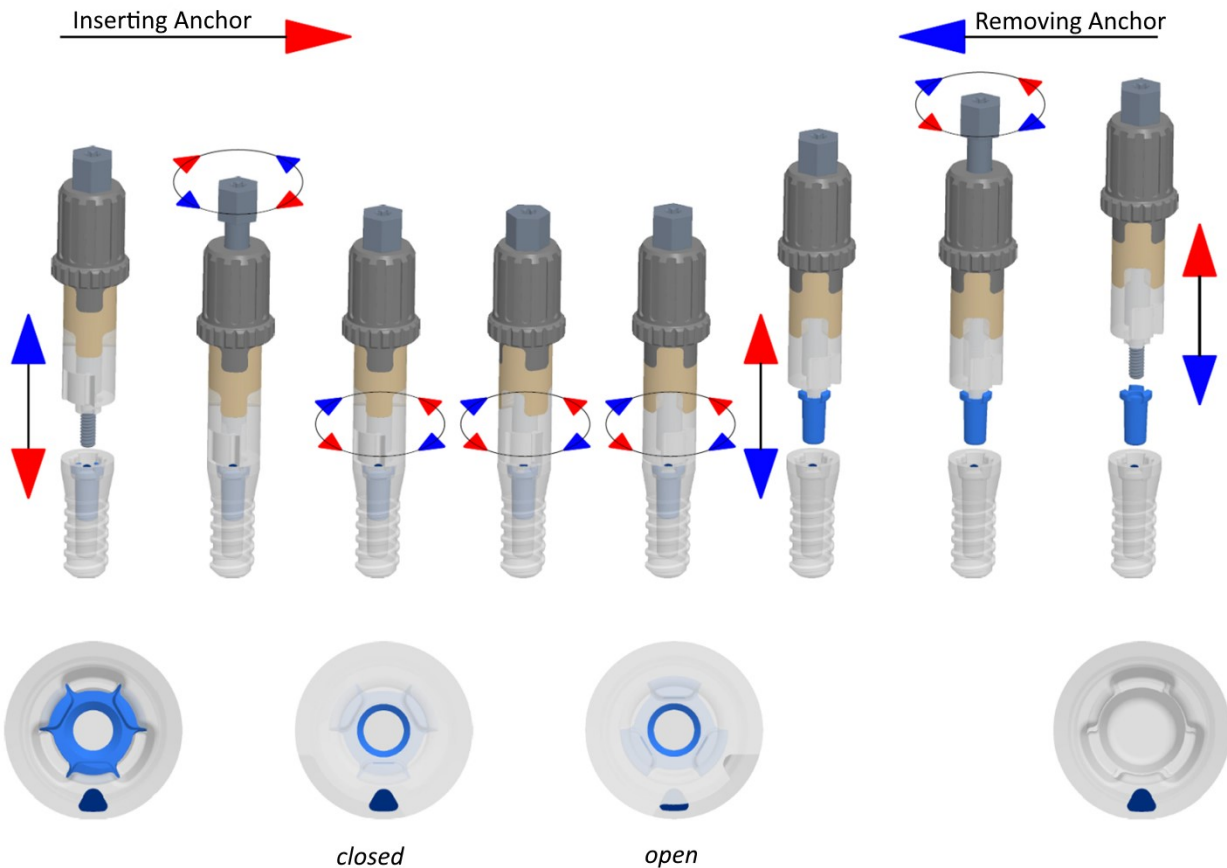
Fig. 2.3-31

Here is an example of how an anchor part can be removed with the instrument components shown. The dot on the implant is right-aligned in the viewing window, and after turning it to the left, the marking moves to the left edge of the viewing window. In this position, the coupled anchor part can be pulled out of the implant.

Anchor part in implant (as supplied)

Only after the impression has been taken and the scaffolding construction has been completed is the pre-assembled anchor part removed from the implant.

With Same Instrument can be an anchor part can be reinserted into the implant if this should become necessary in exceptional cases, e.g. to attach secondary parts.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

2.3.4.7. Planning template

Fig. 2.3-32

Often, an experienced surgeon does not need 3D X-ray diagnostics for surgical planning with the production of drilling templates. If possible, the patient should be spared the considerable radiation exposure caused by CT/CBCT images. In particular, the resolution of dental CBCT images is worse than CT images, and in addition, there is also a transcription error when fixing the drilling templates in the jaw, especially if they are supported by mucous membranes, so that for this reason alone it is not possible to place them with millimeter precision.

In general, however, the planning template we developed for drawing through on printed OPG is very helpful, it is designed for the commonly used magnification factor 1:1.2.

The templates contain a safety area of approx. 2mm on all sides, as distortions and insufficiently accurate magnification ratios must be expected during X-ray diagnostics. In addition, it can be assumed that a height offset error can occur in the case of implantation. In general, an overestimation of the bone supply of up to 6% must be taken into account.

However, it should also be taken into account that, depending on the location, distortions of 30% to 70% in the horizontal direction of some devices can amount to 20% to 30% in the vertical direction, and it must also be taken into account that the magnification factor within an image can fluctuate between a factor of 1.1 and 1.3. In this respect, it is still unavoidable to fix reference bodies of known geometry enorally in the area of desired implant positions and to produce X-ray images, so that a conclusion can be drawn about the distortion factor at the desired implant position from the representation of the reference bodies.

In this context, it is helpful if the planning template is placed over the printed OPG image for each X-ray patient after the implant has been placed, in order to determine how the magnification and distortion factor of the patient's own device is in different locations in order to take these important findings into account for further surgical planning for other patients.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

2.4. Color coding

FARBCODE IMPLANTAT VARIANTEN

Colour Code for Implant Variants								
S	R	W	L	O	OX	O2X	O3X	
Creme	Signal Blue	Jet Black	Heather Violet	Turquoise green	Signal brow	Sulfur yellow	Deep orange	
								
RAL	9001	5005	9005	4003	6016	8002	1016	2011
R	255	0	3	201	15	110	255	237
G	252	46	5	56	112	59	245	92
B	240	122	10	140	51	48	66	41

Signal Red

RAL 3001
R 163
G 23
B 26

Stand: 31.01.2018

ZIRKONUS Implantatsysteme - 73770 Denckendorf - Heerweg 15 d - Tel: +49 (0)711 305 329 - 0 - info@zirkonus.de - www.zirkonus.de 5

Fig. 2.4-01

The color coding makes it easier to assign system components to the implant variants. The implant diameters are color-coded on the implant shoulder.

3. Procedure of the treatment

3.1. Bone preparation

Drilling should always be carried out intermittently (drill the drill should also be pulled completely out of the borehole), whereby drilling should only be carried out for 2 to 3 seconds and under low pressure to the desired depth at speeds according to Fig. 3.1-02 with a maximum drilling depth increase of 1-2 mm. It is important to ensure that the outer area, the drill and the cavity are constantly cooled with pre-cooled (5°C to 8°C) sterile physiological saline solution in order to prevent excessive heating and possible damage to the bone tissue. In addition, bone chips can be removed or rinsed off.

In addition, care must be taken to ensure that the drill does not tilt and get stuck during use (risk of breakage).

At the location of the planned borehole, the drilling site should first be marked with a rose drill by inserting a small trough in order to prevent slipping during the (first) pilot drill.

After the pilot drilling, the core hole diameters are extended to the desired depth between 6 mm and 14 mm in ascending order to the desired core hole diameter between 3 mm and 10 mm with appropriate bone drills and final drills.

The bone bed is then conically shaped with the moulding cutters. This guarantees an optimal fit of the implant in the bone.

A thread is then pre-cut for REGULAR, WIDE and LARGE implants. The SMALL and OVERSIZE implants are self-tapping, no thread needs to be pre-cut.

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

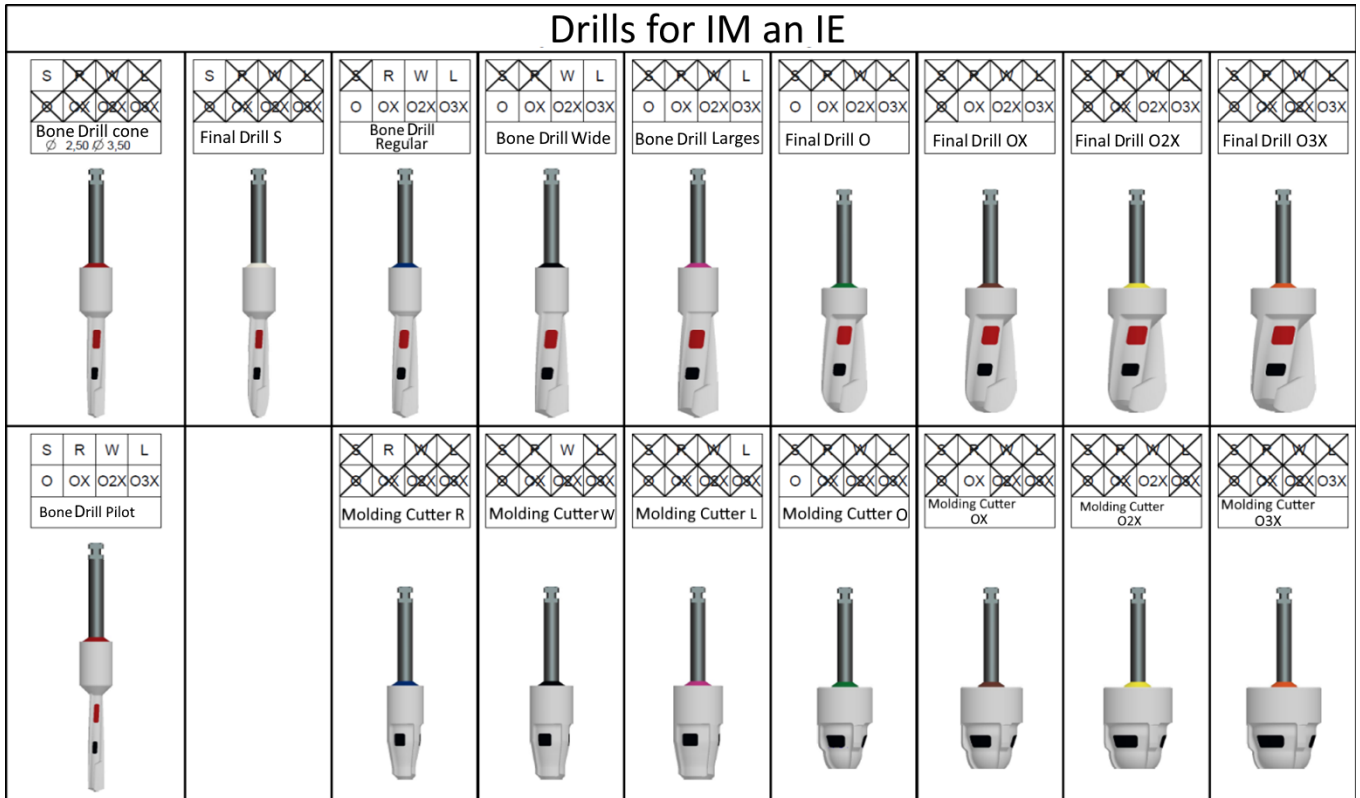
Fig. 3.1-01

Speed recommendations							
	Designation	MatchCode	Speed-recommendation [rpm]		Designation	MatchCode	Speed Recommendation [rpm]
1	Bone Drill Pilot	KB-P	800				
2	Bone Drill Cone	KB-K	800				
3	Final Drill Small	FB-S	600				
4	Bone Drill Regular	KB-R	600		Moulding Cutter R	FF-R	300
5	Bone Drill Wide	KB-W	400		Moulding Cutter W	FF-W	300
6	Bone Drill Large	KB-L	400		Moulding Cutter L	FF-L	250
7	Final Drill O	FB-O	250		Moulding Cutter O	FF-O	250
8	Final Drill OX	FB-OX	200		Moulding Cutter OX	FF-OX	200
9	Final Drill O2X	FB-O2X	200		Moulding Cutter O2X	FF-O2X	200
10	Final Drill O3X	FB-O3X	200		Moulding Cutter O3X	FF-O3X	200

Fig. 3.1-02

See also the instructions for use and preparation attached to the appendix.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

3.2. Detachment from the packaging / attachment of the tool in the case of one-piece and multi-part systems

This is followed by the removal of the sterilized implant body from the blister packaging.

In the one-piece system, the screw-in adapter made of PEEK is pre-mounted on the implant. The ratchet adapter must be attached to the screw-in adapter with light pressure (see Fig. 3.2-01). The implant can then be removed from the implant holder and screwed into the prepared bone cavity.

In the multi-part system, the screw-in tool, consisting of ratchet adapter and screw-in adapter, is placed on the implant fixed in the implant holder (see Fig. 3.2-02) so that the recess of the screw-in adapter lies on the nose of the implant holder and the marking of the implant.

Then the coupling bolt is pushed through into the insertion tool from above and screwed to the anchor part located in the implant until the insertion tool is completely coupled to the implant.

Without having to touch the implant with your fingers, it can be pulled out of the holder sideways with a jerk.

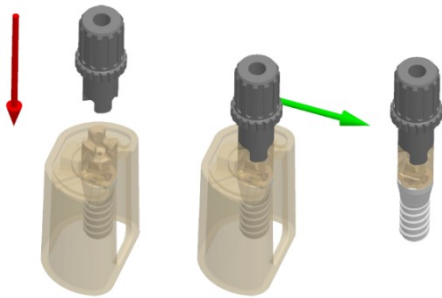
Fig. 3.2-01 One-piece system

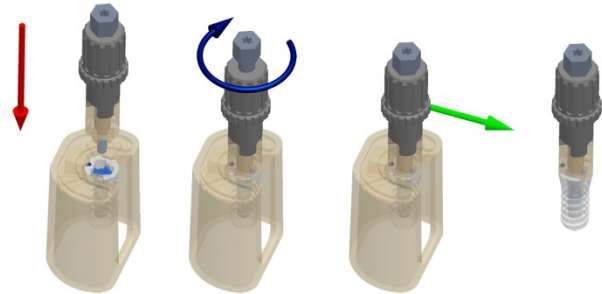
Fig. 3.2-02 Multi-part system

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

3.3. Insertion of the implant into the prepared bone cavity

The implant insertion tool – complex is inserted into the pre-prepared bone storage and then the implant is screwed into the final position using the ratchet. The maximum insertion torque of 35 Ncm must not be exceeded.

The screw-in tool is then decoupled and removed.

Fig. 3.3-01 One-piece system Fig. 3.3-02 Multi-part system

Subsequently, the implant can be closed with a cap for subgingival healing or with a gingival abutment for transgingival healing after insertion.

See Figures 2.3-17, 2.3-18, 2.3-19 in [Chapter 2.3.3.](#)

For wound care, it is important not to use any keratinized mucous membrane – for example by punching the mucous membrane – but to transfer it to the vestibular by means of mucosal flap plasty.

Monofilament or braided suture material is recommended for fixation.

3.4. Exposing

After the healing period/healing phase (usually 6 to 12 weeks), the implant is exposed with the removal of the sealing cap or the healing cap.

See figures in [Chapter 2.3.3.](#)

3.5. Support / fastening of the temporary supply:

The basic planning elements are suitable for attaching the temporary/waxup. The basic planning element is designed in such a way that the mucous membrane is kept away from the implant shoulder.

Until the ceramic abutments / scaffolds are delivered, a temporary/waxup is made for the patient on the planning base elements (PBE). This is the time at which the patient loses his removable temporary, which has been supported on the mucous membrane since then. The new fixed temporary prosthesis is now supported by the implants that have been placed.

See Figure 20 in [Chapter 2.3.3.](#)

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

3.6. Model production / procedure of the production of a superstructure

Insertion of an impression post

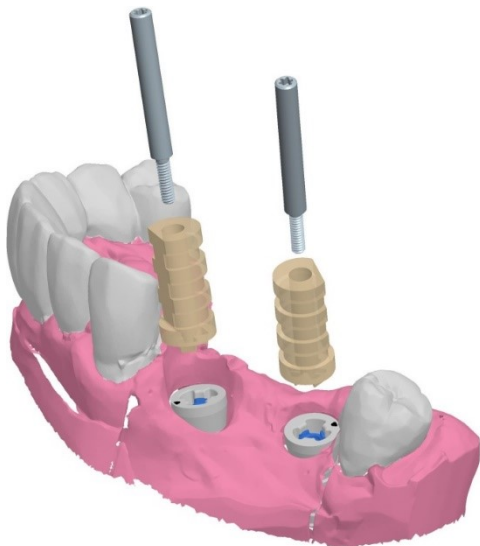
Fig.3.6-01

mounted scanbody

Fig. 3.6-02

After an individual spoon has been produced, an impression of the jaw situation with the implant positions is taken using a conventional impression technique with plastic impression material.

After the plaster model has been made using the ZIRKONUS manipulation implants, the scan bodies are placed and coupled on the manipulation implants located in the plaster model by means of web bolts.

Now the 3D model creation is carried out in the usual way with a suitable scanner (Imetric, ZFX, Evolution or 3Shape 900 series).

First of all, the external geometry of the desired bridge construction must be determined with a suitable dental planning software (e.g. Exocad, etc.). This data must then be transmitted to the company ZIRKONUS Implantatsysteme via ZIRKONUS UPLOADER using the website www.zirkonus.de by data transfer. From there, the data is transferred to the ZIRKONUS planning program and further processed.

The geometry of the individual abutments as well as that of the crowns and bridge construction, including the framework holders, are finally planned at ZIRKONUS and prepared for subsequent production in all-ceramic (zirconium oxide). An intermediate impression is no longer necessary.

Example of a standard abutment:

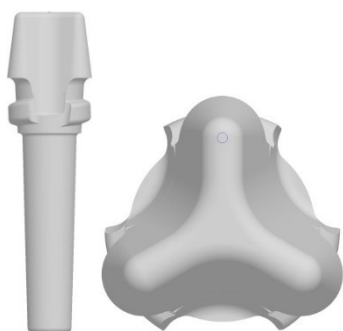
Fig. 3.6-13

Example of an individual abutment:

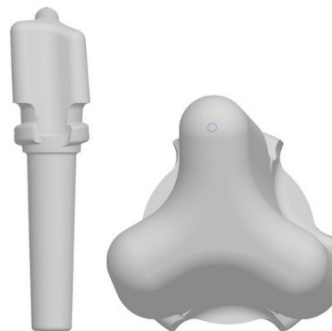
Fig. 3.6-14

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

Abb. 3.6-03 Planung

Fig. 3.6-04 Production

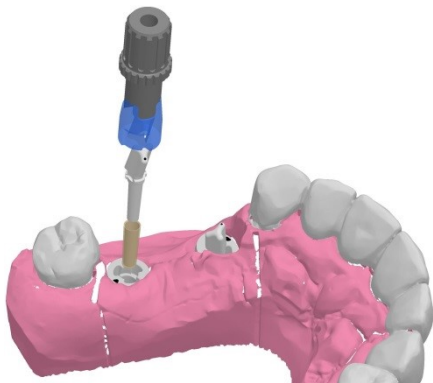
Fig. 3.6-05 Abutment and sleeve

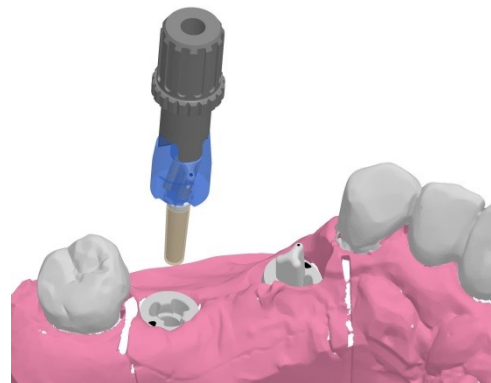
Fig. 3.6-06 Insertion

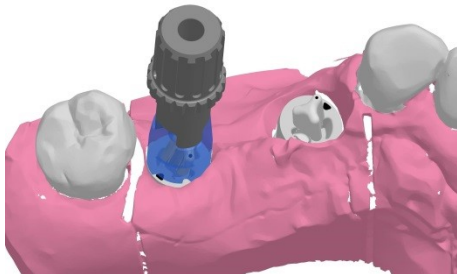
Fig. 3.6-07 inserted, before final position

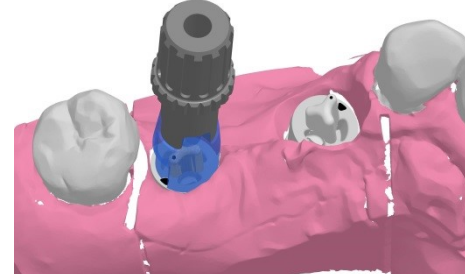
Fig. 3.6-08 inserted, in the end position (60°)

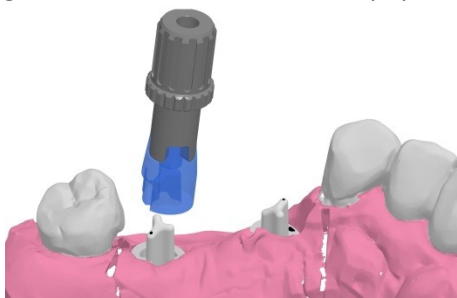
Fig. 3.6-09 Distance adapter

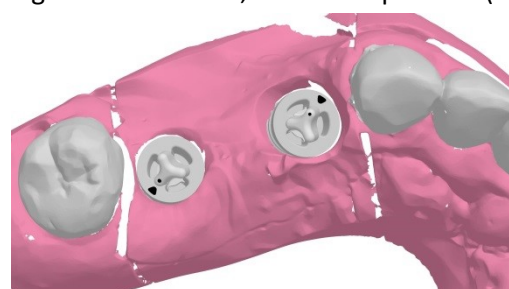
Abb. 3.6-10 Abutment in Endposition

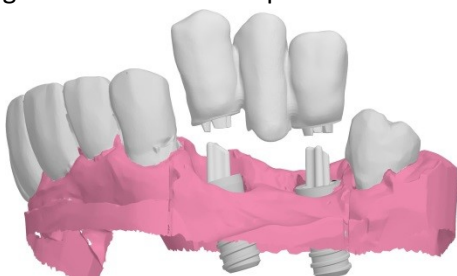
Fig. 3.6-11 Assembly of the superstructure

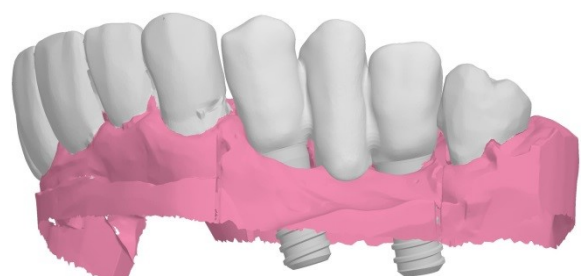
Fig. 3.6-12 Superstructure inserted

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

3.7. Insertion of the abutments / scaffold construction, replacement of the anchor part

After completion of the abutment and scaffolding construction, the anchor parts must first be removed from the implants.

To this end, the following procedure must be followed:

- Removal of the temporary solution
- Removal of the planning base elements from the implants (loosen screw connections)
- Attaching the anchor unscrewing tool and screwing in the coupling bolt for a firm coupling to the anchor part, see Fig. 3.7-01.
- Left rotation by 60° (see marking of the viewing window)
- If the **triangle mark** on the top of the implant shoulder is partially concealed on the left edge of the viewing window in the viewing window, the anchor unscrewing tool can be pulled out with the lever part unlocked, see Fig. 3.7-01.

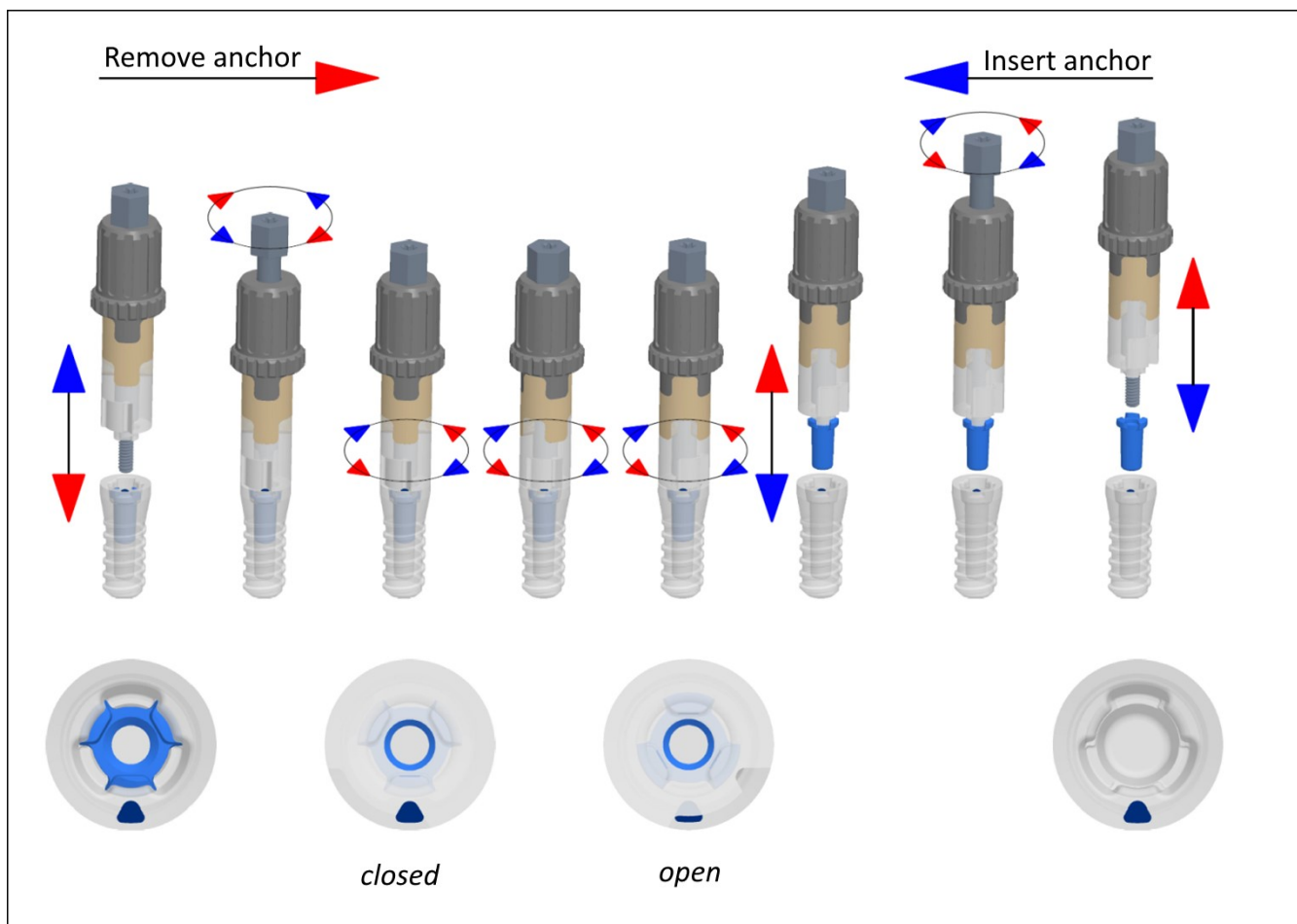
Fig.: 3.7-01

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

3.8. Insertion of the sleeve/abutment/temporary or superstructure

Insertion of the abutment and sleeve using the supplied coupling parts and screw-in tools (please note: the mark must be completely visible in the viewing window (see Fig. 3.8-01)!

Normally, the bridge that has already been made can be inserted without an intermediate impression and the case is thus closed.

In the event that paint or other corrections to the bridge construction are necessary, the following measures will be necessary:

1. Removal of the abutments including sleeves (by means of coupling part, back rotation, left rotation by 60 degrees)
2. Reinsertion of the anchor parts, see Fig. 3.8-01 → Implementation in reverse order!
3. Reattachment of the planning base elements (PBE) or sealing caps in the case of removable interim prosthetics
4. Insertion of the temporary bridge restoration or insertion of the removable interim prosthetics
5. Re-appointment of the patient for definitive insertion of the abutments and bridge construction after correction

(see 3.7. Insertion of the abutments / scaffold construction)

The illustration shown below is the procedure for the multi-part system!

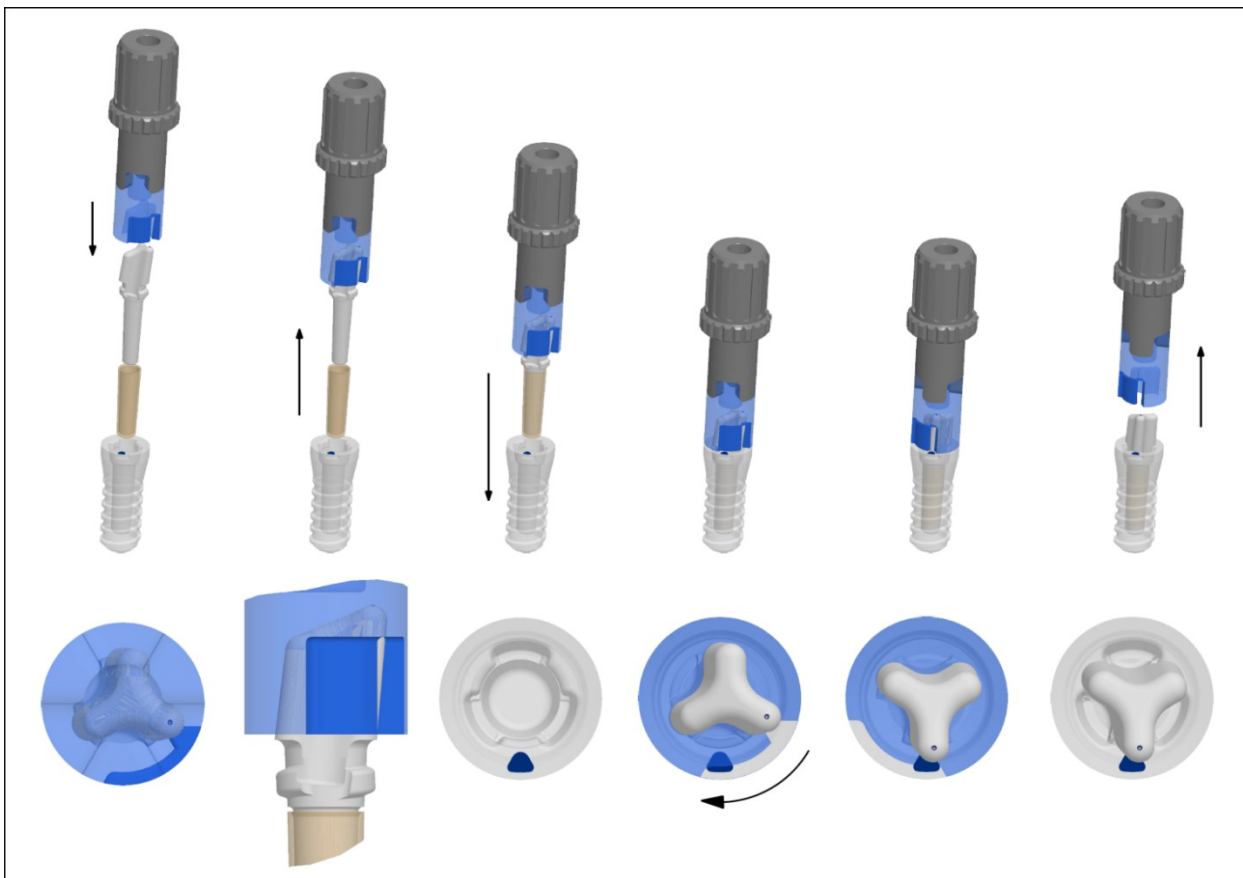
Fig.: 3.8-01

Colour marking fully visible in the viewing window
Abutment in correct position

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

With the **one-piece** system, there is no need to insert the abutment, as the implant body and abutment already form a single unit.

4. Reprocessing of surgical instruments and system components

Surgical instruments and drills supplied by ZIRKONUS are delivered non-sterile and are reusable. The products must be subjected to a complete reprocessing process before their first and each subsequent application.

Prosthetic and system components supplied by ZIRKONUS are also delivered non-sterile and are designed for single use. These products must not be reused and must be subjected to a complete reprocessing process once before use.

For clear identification of the products, please note the information and symbols on the labelling as well as the respective information in the product-specific instructions for use, which is included with each product. For information on validated processes for the preparation of the products, please refer to the corresponding product-specific instructions for use, which are included with each product:

Product group	Delivery form	Again-applicability	Preparation	Directions for use
ZIRKONUS Prosthetic and System Components (Abutments, healing caps, screws, caps, sleeves)	Unsteril 	Single use 	<ul style="list-style-type: none"> Machine cleaning Thermal disinfection Sterilisation (Dampf) 	A 201 - Instructions for use implant system
ZIRKONUS Drill Bit (bone drills, moulding cutters, final drills)	Unsteril 	Max 20 treatment cycles	<ul style="list-style-type: none"> Machine cleaning Thermal disinfection Sterilisation (Dampf) 	A 200 - Instructions for use and preparation
ZIRKONUS surgical instruments (Adapter, Twist Wrench, Thread Cutter)	Unsteril 	Max 20 treatment cycles	<ul style="list-style-type: none"> Machine cleaning Thermal disinfection Sterilisation (Dampf) 	A 200 - Instructions for use and preparation

5. Transport and storage of implants

Transport and storage of the packaged sterile goods is dust, moisture and recontamination protected.

The delivery condition of the ZIRKONUS implants is sterile. The implants are sterilely packaged in a special blister and colour-coded. The selected implant is only removed from the sterile packaging immediately before insertion.

The implants must be stored in their intact sealed packaging in a dry place protected from light and moisture under the following storage conditions:

Instructions for use: A 201 - Instructions for use implant system

- Temperatur: +10°C bis +30°C
- Relative humidity [RH]: 30% - 65%

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

6. Marking / Symbols

Manufacturer within the meaning of Directive 93/42/EEC

Note: NOT STERILE

LOT number

Article

Follow the instructions for use

Attention: Observe warnings!

Note: Sterilized with ethylene oxide

Expiration date

Note: do not reuse!

Note: do not resterilize!

Note: protect from moisture

Note: Protect from direct sunlight

Note: do not use in case of damaged sterile packaging

Note: Storage Conditions Temperature

Note: Storage conditions Humidity

CE mark with identification number of the notified body

0483

mdc medical device certification GmbH, Kriegerstraße 6, 70191 Stuttgart

ZIRKONUS Implantatsysteme GmbH & Co. KG
Bahnhofstraße 18
D-71034 Böblingen / Germany

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

7. Specific error possibilities when dealing with the system

The colour marking on the implant shoulder is a new development. It serves as a guide for the insertion of auxiliary parts and definitive prosthetic parts. With the help of coupling parts, for example, abutments are brought into the correct position. In the viewing window, it is then clearly visible when the correct position of the abutment has been achieved with a 60° rotation, see also Figs. 3.6-05 to Figs. 3.6-10 and also Figs. 3.8-01. It is important that this is done exactly. The use of magnifying glasses is recommended. An instrument with a viewing window display is also available for the anchor part, see Fig. 3.7-01.

Furthermore, it is essential to ensure that while screws are screwed in, the auxiliary part/attachment part lies firmly on the implant shoulder, otherwise there is a risk that the anchor part inside the implant is turned into the wrong position (see Fig. 5.0-01). If the anchor part is in the extension position, it must be rotated clockwise by 60° to the correct position (Fig. 5.0-02) using the anchor ON-OFF rotary tool.

Incorrect seat of the anchor part

(Anchor part unlocked)

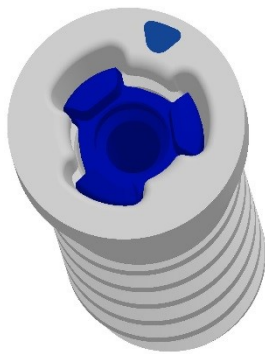
Fig.: 5.0-01

Correct fit of the anchor part

(anchor part locked)

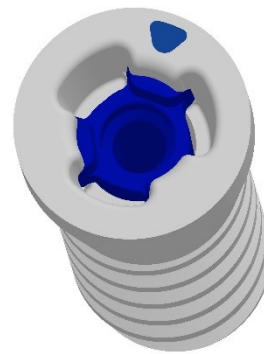
Fig.: 5.0-02

Please note:

- Appropriate measures must be taken to avoid the risk of aspiration when using all system components and prosthetic parts.
 - Never bring ceramic parts together with steel parts, not even when cleaning, disinfecting and sterilizing.
 - No use of force if, for example, parts do not fit together exactly, e.g. because mucous membrane has become trapped between the implant and the gingival former.
 - For tightening and loosening the ceramic screws, e.g. for fastening HDK and PBE, only the PEEK BIT may be used, never the ceramic BIT
 - Ceramic screws can only be screwed in and out with the PEEK bits.
 - For impressions with impression posts, the steel screws "Screw Impression Post" (S-AFP) should be used.
- Never use the BIT rotary wrench in the ceramic version for metal screws! This is only to be used in emergencies to loosen ceramic screws that are stuck due to concretions.
- The metal screw does not match the ZrO2 screw for attaching the HDK & PBE.
 - When taking impressions, when positioning the PEEK impression posts, care must be taken to ensure that they sit on the implant shoulder without gaps, i.e. there must be no mucous membrane trapped

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

between the components. This is the only way to ensure that the 3D model later matches the situation in the jaw. The positioning of the impression post (flattening aligned to the color dot on the implant shoulder) ensures that the exact alignment of the manipulating implant on the model corresponds to the situation in the jaw. This later makes it easier to install the individual abutments in the respective implant.

- We recommend having the zirconium superstructure fully anatomical in the molar area, reduced in the front area for cosmetic reasons for the burning of a glass-ceramic because of the better translucency.
- If the peri-implant bone resorption exceeds 3 mm above shoulder height, the patient must be informed by the practitioner that there is an incalculable risk that the implant construction will break, fall out, or aspirate or be swallowed.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

8. Imprint

Company:

ZIRKONUS Implantatsysteme GmbH & Co. KG

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D-71034 Böblingen

Permanent business establishment:

ZIRKONUS Implantatsysteme GmbH & Co. KG

Heerweg 15D

D-73730 Denkendorf

Phone: +49 (0)711 305329 – 0

Fax: +49 (0)711 305329 – 29

Email: info@zirkonus.de

Website: www.zirkonus.de

District Court Stuttgart HRA 736740

Managing Director: Dr. Dr. Gerd Axel Walther

Top-IdNr.: DE329492060

Legal Notices

No advice

The information provided on the website of ZIRKONUS Implantatsysteme GmbH & Co. KG is not a substitute for individual advice. Our staff will be happy to assist you in person, by phone or by email.

No guarantee for content

The information and information on our websites does not constitute an express or implied assurance or guarantee, in particular not with regard to the nature or suitability for specific purposes.

Trademark and Licensing Rights

Some of the products and designations mentioned on the websites are protected by patents, trademarks and copyrights. The absence of a specific indication or signs [™] does [®] not mean that there is no protection. These websites do not grant third parties any rights to the intellectual property of ZIRKONUS Implantatsysteme GmbH & Co. KG.

Copyrights

The content and structure as well as the texts, images, graphics, files, etc. used on the website are subject to copyright and other laws for the protection of intellectual property. Their disclosure, modification, commercial use or use in other websites or media is not permitted or requires the prior consent of the. All rights reserved.

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9. Grounds:

9.1. Overview Implant Variants and Sizes

		Variants/Sizes															
		S		R		W		L		O		OX		O2X		O3X	
endosseous depth	min. 14 - 15,3 mm	IE-SBE L14	IM-R L14	IE-R L14	IM-WL14	IE-WL14	IM-L L14	IE-L L14	IM-O L14	IE-O L14	IM-OX L14	IE-OX L14	IM-O2X L14	IE-O2X L14	IM-O3X L14	IE-O3X L14	
	min. 11 - 12,3 mm	IE-SBE L11	IM-R L11	IE-R L11	IM-WL11	IE-WL11	IM-L L11	IE-L L11	IM-O L11	IE-O L11	IM-OX L11	IE-OX L11	IM-O2X L11	IE-O2X L11	IM-O3X L11	IE-O3X L11	
endosseous depth	min. 8 - 9,3 mm		IM-R L08	IE-R L08	IM-WL08	IE-WL08	IM-L L08	IE-L L08	IM-O L08	IE-O L08	IM-OX L08	IE-OX L08	IM-O2X L08	IE-O2X L08	IM-O3X L08	IE-O3X L08	
	min. 6 - 7,3 mm				IM-WL06	IE-WL06	IM-L L06	IE-L L06	IM-O L06	IE-O L06	IM-OX L06	IE-OX L06	IM-O2X L06	IE-O2X L06	IM-O3X L06	IE-O3X L06	

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.2. Overview of multi-part implants

ZIRKONUS Implantatsysteme

Implants

Multi-part system									
Size-Variants									
	Ø5,50 Ø4,60 Ø4,05	Ø6,50 Ø5,60 Ø5,05	Ø7,50 Ø6,60 Ø6,05	Ø8,50 Ø7,70 Ø6,90	Ø9,50 Ø8,70 Ø7,90	Ø10,50 Ø9,70 Ø8,90	Ø11,50 Ø10,70 Ø9,90	Plate Ø Cylinder Ø Drill Ø	
	Wide			Large	Oversize	Oversize X	Oversize 2X	Oversize 3X	Min. Enossal
	Regular								
	IM-R L08 (5605502)	IM-W L06 (5606501)	IM-L L06 (5607501)	IM-O L06 (5608501)	IM-OX L06 (5609501)	IM-O2X L06 (5610501)	IM-O3X L06 (5611501)		6 mm
	IM-R L08 (5605502)	IM-W L08 (5606502)	IM-L L08 (5607502)	IM-O L08 (5608502)	IM-OX L08 (5609502)	IM-O2X L08 (5610502)	IM-O3X L08 (5611502)		8 mm
	IM-R L11 (5605503)	IM-W L11 (5606503)	IM-L L11 (5607503)	IM-O L11 (5608503)	IM-OX L11 (5609503)	IM-O2X L11 (5610503)	IM-O3X L11 (5611503)		11 mm
	IM-R L14 (5605504)	IM-W L14 (5606504)							14 mm
									Colour Code
									

This **Surgical Manual** is an **Instruction for use and handling** of the ZIRKONUS implant system

9.3. Overview Implants One-Piece

ZIRKONUS
Implantatsysteme

Implants

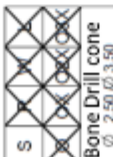
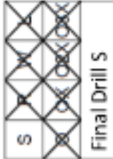
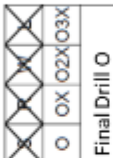
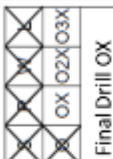
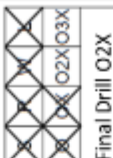
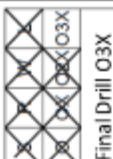
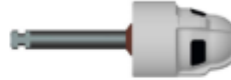
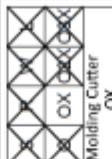
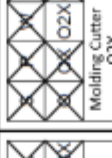
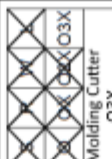
Single-part System										
Size-Variants										
									Plate Ø Cylinder Ø Drill Ø	
Ø3,50	Ø3,50	Ø3,50	Ø5,50	Ø6,50	Ø7,50	Ø8,50	Ø9,50	Ø10,50	Ø11,50	min. Enossal length
Ø3,50	Ø3,50	Ø4,60	Ø5,60	Ø6,60	Ø7,70	Ø8,70	Ø9,70	Ø10,70	Ø11,70	
Ø3,20	Ø3,20	Ø4,05	Ø5,05	Ø6,05	Ø6,90	Ø7,90	Ø8,90	Ø9,90	Ø10,90	
Small Bridge Extension	Small Single Tooth	Regular	Wide	Large	Oversize	Oversize X	Oversize 2X	Oversize 3X		
IE-SBE L11 (5603508)			IE-W L06 (5606506)	IE-L L06 (5607506)	IE-O L06 (5608506)	IE-OX L06 (5609506)	IE-O2X L06 (5610506)	IE-O3X L06 (5611506)		6 mm
IE-SBE L14 (5603509)		IE-R L08 (5605507)	IE-W L08 (5606507)	IE-L L08 (5607507)	IE-O L08 (5608507)	IE-OX L08 (5609507)	IE-O2X L08 (5610507)	IE-O3X L08 (5611507)		8 mm
	IE-SST L11 (5603510)	IE-R L11 (5605508)	IE-W L11 (5606508)	IE-L L11 (5607508)	IE-O L11 (5608508)	IE-OX L11 (5609508)	IE-O2X L11 (5610508)	IE-O3X L11 (5611508)		11 mm
	IE-SST L14 (5603511)									14 mm
										Colour Code

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.4. Overview of rotating instruments

Drills

ZIRKONUS
Implantatsysteme

 <p>Bone Drill cone ∅ 2,50 ∅ 3,50</p>		 <p>Bone Drill Pilot</p>		 <p>Final Drill S</p>		 <p>Bone Drill Regular</p>		 <p>Bone Drill Wide</p>		 <p>Bone Drill Large</p>		 <p>Final Drill O</p>		 <p>Final Drill OX</p>		 <p>Final Drill O2X</p>		 <p>Final Drill O3X</p>		 <p>Molding Cutter O</p>		 <p>Molding Cutter OX</p>		 <p>Molding Cutter O2X</p>		 <p>Molding Cutter O3X</p>
--	---	--	---	--	---	---	---	--	---	---	--	---	---	---	---	--	---	--	---	--	---	---	---	--	---	--

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system


9.5. Overview of depth gauges and thread cutters

ZIRKONUS Implantatsysteme

Gauges and Threadcutter

Gauges								
	TL-S	TL-R	TL-W	TL-L	TL-O	TL-OX	TL-O2X	TL-O3X
40 mm								

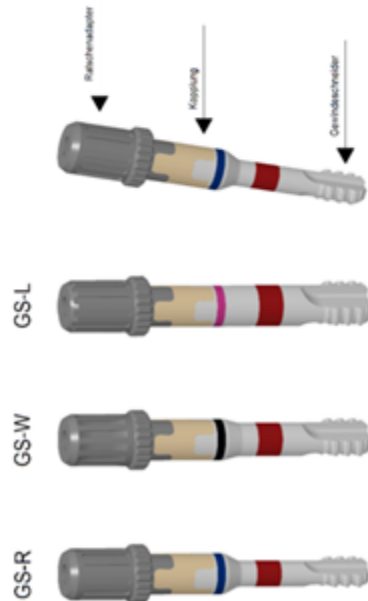Threadcutter

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.6. Overview Cap - Gingival Abutment

Cap/Gingiva Molder

	S	R	W	L	O	OX	O2X	O3X
Cap HDK IM V								
Collared Gingiva Molder HDK IM G								
Zylindrical Gin- giva molder 2,5 mm HDK IM Z025								
Zylindrical Gin- giva molder 4,50 mm HDK IM Z045								
IE-System (PEEK)	SBE SST							

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.7. Overview of planning basic elements

Base-Elements

IM-System (PEEK)	Base Element PBE-IM-	S	R	W	L	O	OX	O2X	O3X
IE-System (PEEK)	Base Element PBE-IE-	SBE SST							

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.8. Overview of impression posts

ZIRKONUS

Implantatsysteme

Molding Post (AFP)

		Variants/Sizes							
		S	R	W	L	O	OX	O2X	O3X
IM-System			AFP-IM-R	AFP-IM-W	AFP-IM-L	AFP-IM-O	AFP-IM-OX	AFP-IM-O2X	AFP-IM-O3X
IE-System	SBE	AFP-IE-SBE	AFP-IE-R	AFP-IE-W	AFP-IE-L	AFP-IE-O	AFP-IE-OX	AFP-IE-O2X	AFP-IE-O3X
	SST	AFP-IE-SST							

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.9. Overview of impression posts on implants

ZIRKONUS

Implantatsysteme

Molding Post on Implant

		Variants/Sizes							
		S	R	W	L	O	OX	O2X	O3X
IM-System			AFP-IM-R	AFP-IM-W	AFP-IM-L	AFP-IM-O	AFP-IM-OX	AFP-IM-O2X	AFP-IM-O3X
IE-System	SBE	AFP-IE-SBE	AFP-IE-R	AFP-IE-W	AFP-IE-L	AFP-IE-O	AFP-IE-OX	AFP-IE-O2X	AFP-IE-O3X
	SST	AFP-IE-SST							

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.10. Overview of Manipulating Implants

Rigging-Implants


		Variants/Sizes							
		S	R	W	L	O	OX	O2X	O3X
IM-System			M-IM-R	M-IM-W	M-IM-L	M-IM-O	M-IM-OX	M-IM-O2X	M-IM-O3X
IE-System	M-IE-SBE		M-IE-R	M-IE-W	M-IE-L	M-IE-O	M-IE-OX	M-IE-O2X	M-IE-O3X
	M-IE-SST								

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.11. Insertion process

ZIRKONUS Implantatsysteme

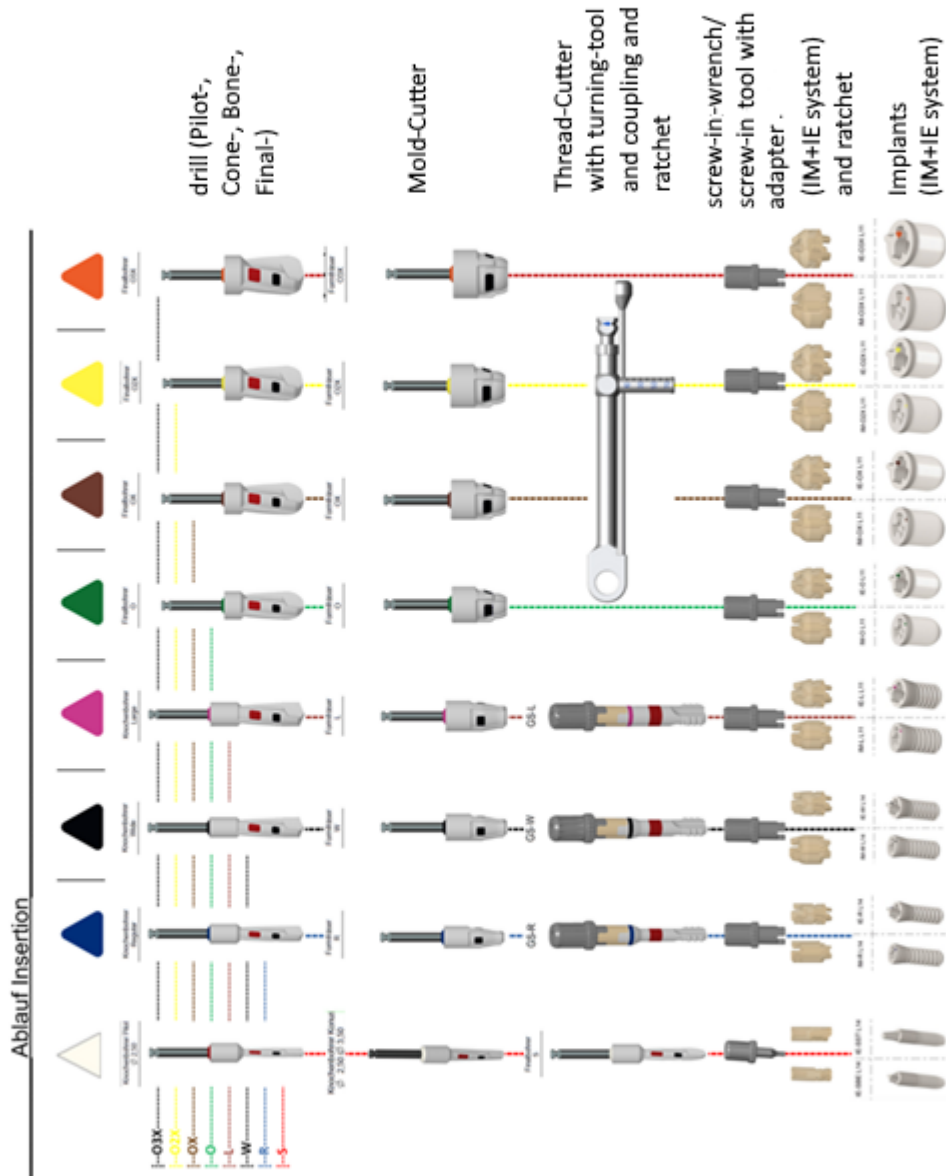Sequence of Insertion

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.12. Overview Scanbody

ZIRKONUS

Implantatsysteme

Scanbody

		Variants/Sizes							
		S	R	W	L	O	OX	O2X	O3X
IM-System			SB-IM-R	SB-IM-W	SB-IM-L	SB-IM-O	SB-IM-OX	SB-IM-O2X	SB-IM-O3X
IE-System	SB-IE-SBE		SB-IE-R	SB-IE-W	SB-IE-L	SB-IE-O	SB-IE-OX	SB-IE-O2X	SB-IE-O3X
	SB-IE-SST								

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.13. Overview Scanbody on Manipulating Implant

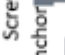
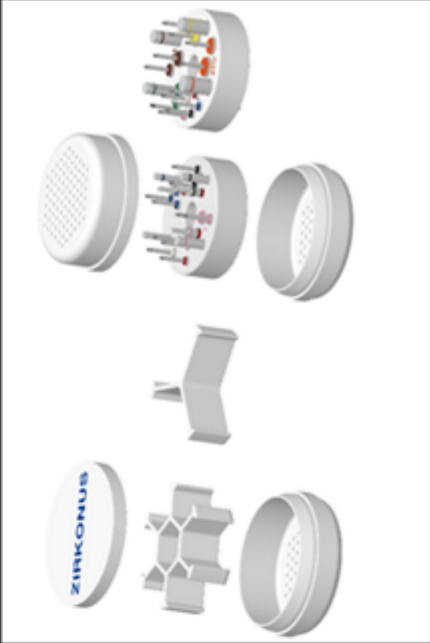
Scanbody on Rigging-Implant

		Variants/Sizes							
		S	R	W	L	O	OX	O2X	O3X
IM-System			SB-IM-R	SB-IM-W	SB-IM-L	SB-IM-O	SB-IM-OX	SB-IM-O2X	SB-IM-O3X
IE-System		SB-IE-SBE	SB-IE-R	SB-IE-W	SB-IE-L	SB-IE-O	SB-IE-OX	SB-IE-O2X	SB-IE-O3X
		SB-IE-SST							

This *Surgical Manual* is an *Instruction for use and handling* of the ZIRKONUS implant system

9.14. Overview of System Components

System components

		Anchor	Shell	Abutment S	Abutment I	Screw	Screw	Healing phase	Final Situation	Template for planning
Parts, short	Name	 ANKER-K	 HUE-K	 ABUT-STAND-K	 ABUT-INDIV-K	 S-Z16-L075 Innendrehrund	 S-Z16-L085 Innendrehrund	 Implant, Anchor, Cap, Ceramic-screw, multi-part system	 Implant, Shell, Abutment, Crown multi-part system	 ZIRKONUS Implantatsysteme
Parts, long	Name	 HUE-L	 ABUT-STAND-L	 ABUT-INDIV-L	 S-T6-L085 Innensechsrund	 single-part system	 single-part system	 Primary packaging (Blister)	
Wrench	 MA-SW4	 Screw Anchor	 screw-in-tool	 screw-in wrench	 ZIRKONUS	 ZIRKONUS	 Adapter turning-key with RA-shaft			

This Surgical Manual is an Instruction for use and handling of the ZIRKONUS implant system

9.15. A 200 - Instructions for use and preparation

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